



ORMESA®
• MADE IN ITALY SINCE 1980 •

DONDOLINO



DONDOLINO

Stability in motion.

Dondolino is a prone standing frame designed and built by Ormesa to help children acquire **independence in maintaining a standing posture**.

Dondolino's use in everyday contexts **encourages environmental interaction** by facilitating upper limb and head mobility and the development of individual motor skills.

Dondolino is an extremely versatile tool that comes in three sizes with a wide range of adjustments. It can be used by children and teens from 75 to 150 cm tall, **adapting to their growth, development of individual skills and hip and knee movements**.



PRODUCT
REGISTERED



PRODUCT
MARKED







Tilting frame

Through spring stopper adjuster it is possible to switch from the upright position to a more prone one in order to stimulate the active control of the head and trunk.

These changes of attitude can be applied in a completely fluid way even with the child on board, without causing pathological reactions of postural discomfort.

Dondolino contains without forcing, adapting the structure to different postural needs.



**Technical features
and adjustments**



Frame

for children needing to improve straightening reactions. Protrusion free to guarantee safety. Allows trunk rotation and maximum upper limb freedom of movement. Smooth continuous adjustment (through spring stopper adjuster).



Tools free

adjustment without tools but with special knobs, designed by Ormesa, with effective and safe grip. The knobs effectively maintain tightening through an internal mechanism that does not mark and ruin the tube.




Graduated bars




to help finding the right postural setting chosen by the therapist.



Pelvic support and thoracic support

independently adjustable in circumference, height and forward-backward. If the user has sufficient trunk control, the thoracic support can be removed.

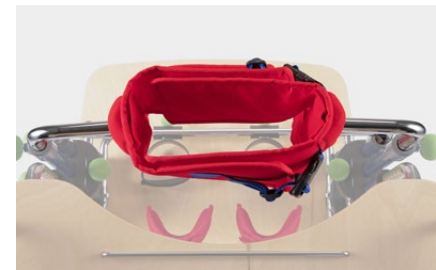

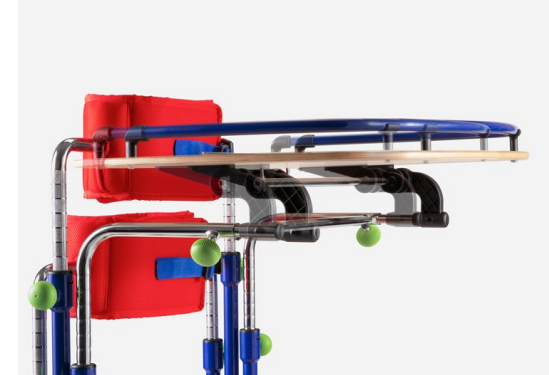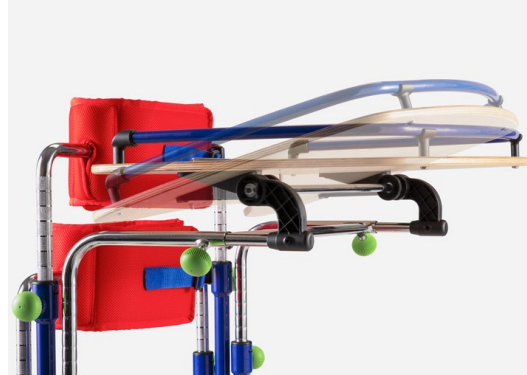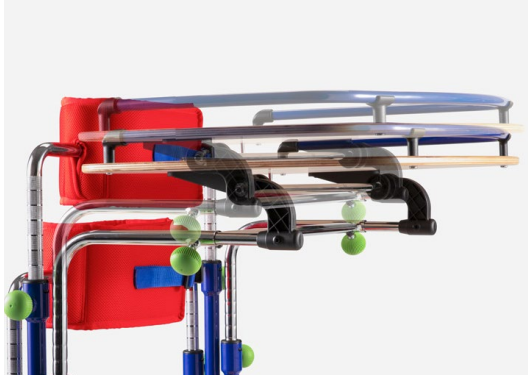

Birch plywood table

designed for educational purposes and for supporting upper limbs, can be removed to facilitate upper limb and trunk exercises. Height, inclination and forward-backward adjustment.




Knee supports

with removable, washable padding and with holes for patella discharge. Independent height, forward-backward, rotation and abductive adjustments.



Heel stops

width and depth adjustments.

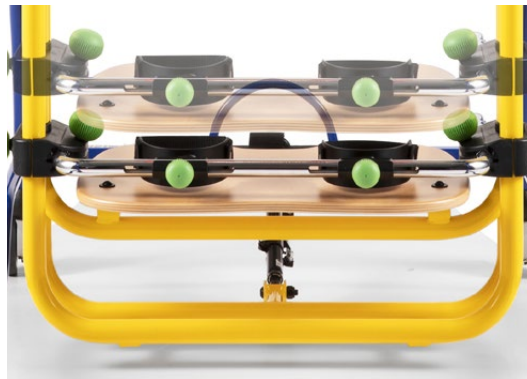

Birch plywood platform

height adjustments.




Twin swivel wheels

75 mm diameter, all with brakes. The wheels allow movement from one room to another, even during use.




+	-	×	÷	=
A	B	C	D	E
F	G	H	I	J
K	L	M	N	O
P	Q	R	S	T
U	V	W	X	Y
Z				

other components

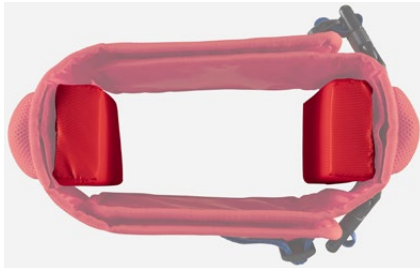

944 Circumference reducers

each 3 cm thick, can be inserted into supports for thinner users.

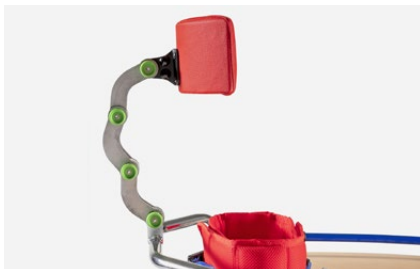

865 Multiadjustable headrest


with height, inclination and forward-backward adjustment knobs.





Dondolino sizes and weight


	size 1	size 2	size 3
A: base encumbrance	63 cm	77 cm	83 cm
B: length encumbrance	from 75 to 85 cm	from 81 to 91 cm	from 85 to 95 cm
C: platform to table height	from 45 to 79 cm	from 52 to 90 cm	from 62 to 107 cm
D: overall width	63 cm	67 cm	71 cm
E: platform to knee pads height	from 12 to 26 cm	from 12 to 36 cm	from 14 to 48 cm
F: platform to pelvic support height	from 39 to 59 cm	from 47 to 74 cm	from 58 to 93 cm
G: platform to thoracic support height	from 50 to 76 cm	from 56 to 90 cm	from 70 to 110 cm
L: circumference of supports	60 cm	65 cm	70 cm
M: user height	from 75 to 100 cm	from 90 to 125 cm	from 120 to 150 cm
P: frame inclination	from 0° to 19°	from 0° to 14°	from 0° to 11°
Total weight	21 kg	24 kg	26 kg
Maximum load	35 kg	45 kg	55 kg

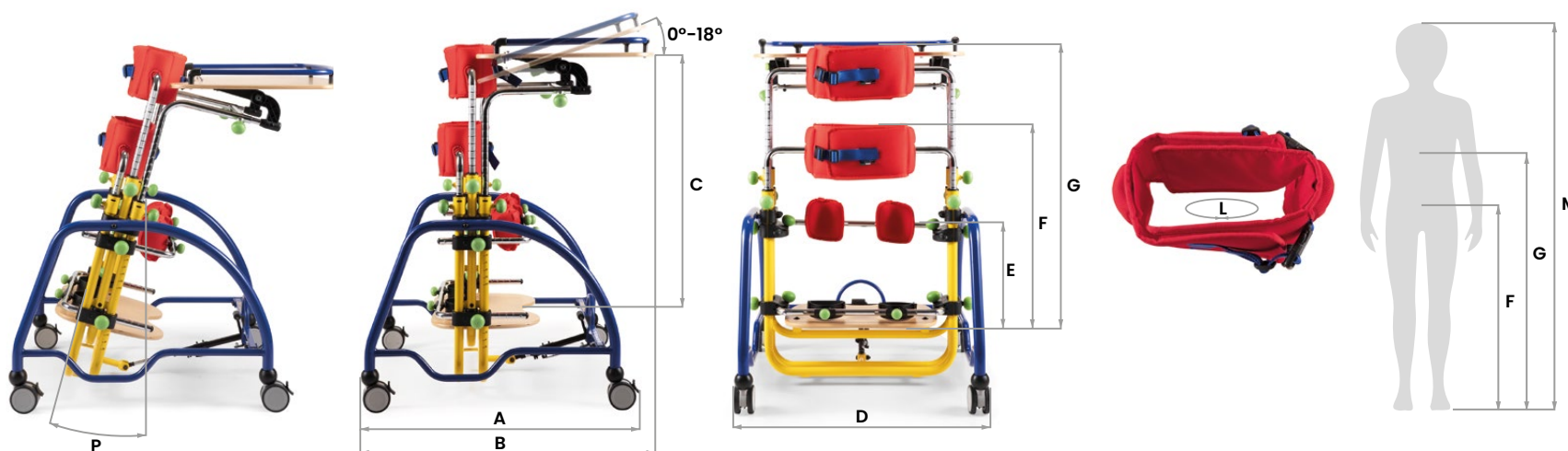

Make sure that measures F and G indicated in the figure do not exceed the suggested range indicated for the supports.

Pietro and his story await you in
official website Ormesa, together
with other protagonists of our
fantastic global community.



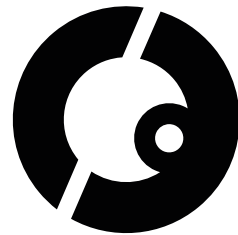

Manufactured by Ormesa srl

Available from Medifab Ltd in New Zealand.
Call 0800 543 343 | solutions@medifab.com

www.medifab.co.nz