

# Houdini 13 Harness

- Can be fitted on **wheelchairs and some easy chairs** as well as using in vehicles
- All buckles are **behind the seat** so the user cannot release them



Model	Standard Push Press Buckle	Pen-Key Release Buckle	Age Group
13J	4604-1302-017	4604-1302-031	3 to 8 years
13H	4604-1303-017	4604-1303-031	7yrs to Adult

## MODEL 13 BENEFITS

Suitable for use with people who are confused or demonstrate behaviours of concern, as all the buckles are behind the seat back.


In vehicles, this harness is suitable for use on a **single seat only**, and requires a gap between the vehicle seat and back support for installation.

Not a safety harness: must be worn in conjunction with vehicle safety belt.

# Houdini 22 Harness

TAXI

- Can be used on **single or bench seats**
- Fit to the user **before entering vehicle**
- Good option for **taxi's** and also suitable for hip spicas

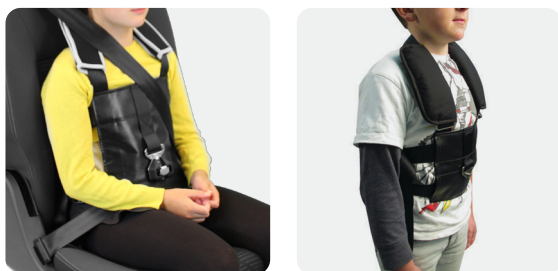


Model	Anti-Escape Buckle	Pen-Key Release Buckle	Age Group
22C	4604-2211-017	4604-2211-031	3 to 8 years
22D	4604-2212-017	4604-2212-031	7yrs to Adult

## MODEL 22 BENEFITS

Consists of a removable torso section which is attached to the individual prior to entering the vehicle.

Once the individual is seated, two buckles are engaged to provide safe travel.

Not a safety harness: must be worn in conjunction with vehicle safety belt.

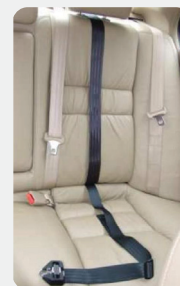

**Additional seat attachment strap for use in another vehicle:**  
4604-2222-017

# Houdini 27 Harness



- **Centre-release buckle** to provide postural support, **anti-escape** to prevent removal
- Attachment straps **fasten horizontally** around seat
- Suitable for single back support seats and **air travel\***



Model	Centre-Release Buckle	Anti-Escape Buckle	Pen-Key Release Buckle	Age Group
27A	4604-2702-017	4604-2702-026	4604-2702-031	3 to 8yrs
27B	4604-2703-017	4604-2703-026	4604-2703-031	7yrs to Adult

## MODEL 27 BENEFITS

Designed to offer postural support and a certain degree of restraint as it is a full 5-point harness. Identical to Houdini 28; however straps fasten horizontally (28's fasten vertically).

Only suitable for single-backed seats.

Not a safety harness: must be worn in conjunction with vehicle safety belt.

\* Signed approval must be obtained from the airline prior to travel.

# Houdini 28 Harness

- Medifab's **most popular** model of Houdini Harness
- Attachment straps **fasten vertically** around seat
- **Centre-release buckle to provide** postural support, **anti-escape** to prevent removal



Model	Anti-Escape Buckle	Centre-Release Buckle	Pen-Key Release Buckle	Age Group
28	-	4604-2803-017	4604-2803-031	Junior to Adult
28SB	4604-2803-026	-	-	Junior to Adult

## MODEL 28 BENEFITS

Designed to offer postural support and a certain degree of restraint as it is a full 5-point harness. This harness is the same as the Model 27, however the seat attachment straps go vertically around the vehicle seat.

This is suitable for **single or bench seats**. For seats that 'narrow' towards the top of the back support, the Model 27 is more suitable. Not a safety harness: must be worn in conjunction with vehicle safety belt.

# Houdini 31 Harness

- **Most secure harness**, with twin anti-escape buckles
- Can be used on **deep 'bucket' single or bench seat**



Model	Anti-Escape Buckle	Pen-Key Release Buckle	Age Group
31A	4604-3101-017	4604-3101-031	3 to 8yrs
31B	4604-3102-017	4604-3102-031	7yrs to Adult

## MODEL 31 BENEFITS

This harness offers more support due to a front PVC pad that is attached around the individual's torso with two straps fitted with anti-escape buckles.

Not a safety harness: must be worn in conjunction with vehicle safety belt.

## CROSS STRAP

Fit to a harness to hold shoulder straps together & prevent them slipping off.

**PVC with side-release buckle**  
4604-3112-019

**PVC with anti-escape**  
4604-3004-026

**PVC with pen release**  
4604-3004-025




**Cross Strap Padded Sleeve (side release buckle)**  
Sml 4604-3117-017  
Lge 4604-3118-017



Designs and styles may vary.

## BUCKLE OPTIONS:

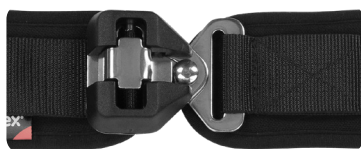

### Anti-Escape

Successfully prevents those with inquisitive fingers, challenging behaviour or learning difficulties from releasing the buckle. Operates in reverse to normal, as it is disengaged by raising a spring-loaded lever and will not press release.



### Pen/Key Release

Requires a tool such as a pointed pen or key to prod the buckle release mechanism to open. Medifab supply a release tool with each of these buckles.