



Snooooooze

This flexible sleep system is designed for children and adults. It is lightweight, extremely portable and highly adjustable to cater for changing positioning needs and growth.

Materials

- The receptor sheet for Snooooooze is made from a soft and comfortable brushed nylon.
- Temperature control fitted sheets are also available which feature a soft terry towelling contact surface, paired with a cooling Outlast underlay for active cooling.
- The supports are filled with a chipped visco-elastic memory foam which provides support and immersion.
- The majority of supports are supplied as standard with two layers of cover; an inner cover made of a fire resistant material, and a standard outer cover made from wipe-down, PU-coated fabric which can be removed and washed.

- The inner covers acts to contain the filling when the outer cover is being washed. This inner cover also features a zip opening, which allows foam to be removed or added, in order to adjust the support volume depending on the user's needs.

Cover Options

- The majority of supports are available with an optional outer cover which can also be removed and washed. This cover is made from a two-layer fabric consisting of a soft and breathable cotton contact surface, with an underlay made from Outlast material which actively cools and assists with lowering temperature.



Support in 3 lying positions

The Snooooooze system supports the user in supine (on their back), side and semi side-lying positions. This key feature enables us to accommodate and correctly position posture overnight, whilst maintaining a high level of comfort.

Supine lying

To ensure correct positioning both during the day and night, we stabilise the pelvis. The pelvis is the foundation of good posture and is therefore, the first area we evaluate when conducting an assessment. Using a component like the W Leg Trough, we aim to draw the legs into mid-line.

Side-lying

The Body Roll has heralded a new stage in side-lying, and can help achieve an excellent trunk and leg position for both left and right side-lying. It is most effective when the trunk has

shortening sides and the legs are windswept. Alignment of the legs is achieved when a pillow is placed between the user's legs. The final key to achieving good posture is to support and position the head.

Sleep for everyone

When the user struggles to sleep, it has a direct impact on the sleep patterns of every family member. Caregivers can be left exhausted, unable to think clearly and struggling to cope with their daily activities. By ensuring that the user has a good night's sleep, it also means that their family will enjoy the many benefits of sleep too.

Receptor Sheet

- For secure positioning, this soft loop mattress sheet ensures that laterals and brackets stay in place, provided they have Velcro on their base. It is highly effective with both Curved Wedges and Tri-Fold Laterals.
- Sizes to suit a variety of mattresses.
- Machine washable.



| | Cot | Single | King Single | Double |
|---------------|---------------|---------------|---------------|---------------|
| Code | 7004-0021-000 | 7004-0022-000 | 7004-0023-000 | 7004-0024-000 |
| Length | 1300mm | 1880mm | 2030mm | 1880mm |
| Width | 690mm | 920mm | 1070mm | 1380mm |

Fitted Temperature Regulation Sheets

- Fitted sheet providing a comfortable and cool layer for sleep.
- Available in a range of mattress sizes.

| | Cot | Single | King Single | Double |
|---------------|---------------|---------------|---------------|---------------|
| Code | 7004-0031-000 | 7004-0032-000 | 7004-0033-000 | 7004-0034-000 |
| Length | 1450mm | 2000mm | 2100mm | 2200mm |
| Width | 800mm | 1000mm | 1080mm | 1500mm |



Tri-Fold Positioning Brackets

- Can be used alone to directly support the body, or with other components from the Snooooooze range, such as the Body Roll.
- Has an internal bracket to ensure shape and support is maintained.
- Removable and washable cover.

| | Paediatric | Small | Large |
|--------------------|---------------|---------------|---------------|
| Code | 7004-0111-000 | 7004-0112-000 | 7004-0114-000 |
| Base Length | 100mm | 150mm | 200mm |
| Height | 85mm | 130mm | 200mm |
| Width | 85mm | 130mm | 200mm |



Curved Wedge

- These curved wedges can be used to directly support the body, or to support other components in the sleep system.
- Crafted to conform to the body's contours, they aid in proper body positioning.
- Velcro attachment secures the wedges to the receptor sheet.

| | Paediatric | Small | Large |
|-----------------------|---------------|---------------|---------------|
| Code | 7004-0131-000 | 7004-0132-000 | 7004-0134-000 |
| Length at Base | 120mm | 125mm | 250mm |
| Height | 100mm | 125mm | 150mm |
| Width at Base | 100mm | 125mm | 150mm |



Banana Cushion

- A multi-functional and versatile piece of sleep equipment. For example, it can be used as a short body roll, under a W Leg Trough, or under the neck.
- Removable and washable cover.
- A thermal regulation cover is available for this product.

| | Paediatric | Large |
|-------------------|-----------------|-----------------|
| Code | 7004-0301-000 | 7004-0304-000 |
| Diameter | 100mm | 180mm |
| Length | 450mm | 850mm |
| Width | 120mm | 200mm |
| Cover Code | 7004-0301-000-C | 7004-0304-000-C |



Horseshoe

- Versatile support which can be used for multiple applications in supine and side-lying.
- A thermal regulation cover is available for this product.
- Removable and washable cover.


| | |
|-------------------|-----------------|
| Code | 7004-0313-000 |
| Length | 1000mm |
| Diameter | 210mm |
| Cover Code | 7004-0313-000-C |



Star Head Support

- Adjustable head support with profiling points for added stability or alignment correction.
- For people with low to medium tone, the Star Head Support is ideal, as each of the points can be folded to create a shape which can be used to position the head as needed.
- Can be supplemented with wedges.
- Removable and washable covers.
- A thermal regulation cover is available for this product.

| | Paediatric | Large |
|-------------------------------|-----------------|-----------------|
| Code | 7004-0801-000 | 7004-0804-000 |
| Length Star Tip to Tip | 560mm | 690mm |
| Diameter Inner Circle | 80mm | 100mm |
| Cover Code | 7004-0801-000-C | 7004-0804-000-C |

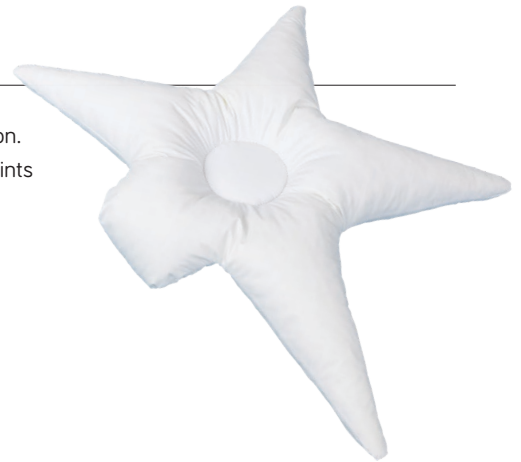

Roll

- The Long Roll cushion can be used as a head, knee or back support.
- The Long Roll is a multi-functional cushion that can be used as a short body roll, under a W leg trough, under a neck or support other areas around the body.
- Removable and washable cover.
- A thermal regulation cover is available for this product.

| | |
|-------------------|-----------------|
| Code | 7004-0322-000 |
| Diameter | 140mm |
| Length | 600mm |
| Cover Code | 7004-0322-000-C |



Body Roll

- Long positioning piece offering full length support for users.
- Beneficial for side and semi side-lying positions, offering support either in front of or behind the user.
- Can also provide support between the legs for protection of the user's top hip joint.
- A thermal regulation cover is available for this product.

| | Paediatric | Large |
|-------------------|-----------------|-----------------|
| Code | 7004-0331-000 | 7004-0334-000 |
| Length | 1300mm | 2000mm |
| Diameter | 90mm | 150mm |
| Cover Code | 7004-0331-000-C | 7004-0334-000-C |



T-Roll

- T-Roll leg support and memory foam cushion supports the user's legs in abduction and inhibits adduction.
- For users who find the W Leg Trough unsuitable due to contracture at the knee, T-Rolls are an ideal solution for a presentation of legs that come to mid-line.
- Can be used with other components of the sleep system range.
- Washable and removable cover.
- A thermal regulation cover is available for this product.

| | Paediatric | Small | Large |
|----------------------------------|-----------------|-----------------|-----------------|
| Code | 7004-0201-000 | 7004-0202-000 | 7004-0204-000 |
| Diameter of Main Roll | 150mm | 200mm | 250mm |
| Diameter of Centre Roll | 150mm | 200mm | 250mm |
| Width Either Side of Roll | 225mm | 250mm | 250mm |
| Height of Centre Section | 150mm | 200mm | 250mm |
| Cover Code | 7004-0201-000-C | 7004-0202-000-C | 7004-0204-000-C |



Curved T-Roll with Wrap

- Same great benefits of the standard T-roll, but with a curved design for improved contouring around the underside of the legs.
- Supplied standard with wrap.
- Washable and removable cover.
- A thermal regulation cover is available for this product.

| | |
|----------------------------------|-----------------|
| Code | 7004-0230-000 |
| Diameter of Main Roll | 250mm |
| Diameter of Centre Roll | 250mm |
| Width Either Side of Roll | 250mm |
| Cover Code | 7004-0230-000-C |



T-Roll with Stabiliser

- Same great benefits of the standard T-roll, but with an added section to improve stability.
- Washable and removable cover.
- A thermal regulation cover is available for this product.

| | Paediatric | Large |
|----------------------------------|-----------------|-----------------|
| Code | 7004-0211-010 | 7004-0214-010 |
| Diameter of Main Roll | 150mm | 250mm |
| Diameter of Centre Roll | 150mm | 250mm |
| Width Either Side of Roll | 225mm | 250mm |
| Cover Code | 7004-0211-000-C | 7004-0214-000-C |



Twin Side Rolls

- The Twin Side Rolls and flap are a positioning element to centralise posture, often used by clients who have a feeding tube and require elevation during feeding.
- Can be used with other components of the sleep system range.
- Washable and removable cover.

| | |
|--------------------------------------|---------------|
| Code | 7004-0820-000 |
| Length of Roll(s) | 900mm |
| Diameter of Roll(s) | 150mm |
| Width at Top Between Rolls | 460mm |
| Width at Bottom Between Rolls | 360mm |



Wedge

- Versatile support wedge for either the head or trunk.
- Where the neck tone is too great for the Star, a Wedge can be used up against the head, usually with the clients own pillow padding. Opposite positioning wedges may be needed on the trunk.
- Can also be used as a support against another component of the sleep system range.
- Removable and washable cover.

| | Paediatric | Small | Large |
|--------------------|---------------|---------------|---------------|
| Code | 7004-0121-000 | 7004-0122-000 | 7004-0124-000 |
| Base Length | 100mm | 200mm | 405mm |
| Height | 80mm | 150mm | 150mm |
| Width | 80mm | 150mm | 150mm |



W Leg Trough With Wrap

- W Leg Trough support cushion positions the legs, preventing rotation and reducing the risk of hip dislocation.
- Leg alignment is achieved with the whole leg supported and positioned correctly. The wrap (can be used with or without the wrap) prevents the client from crossing their legs and holds the position securely.
- **W Leg Trough With Wrap - 3 Cells** - The same great benefits of the W Leg Trough however this version has the added flexibility to adjust the filling within each area of the cushion independently, tailoring the support for the user. It provides opportunity to adjust for positions such as abduction, adduction, rotation and windsweeping.
- Removable and washable covers.
- A thermal regulation cover is available for this product.

| | Paediatric | Small | Medium | Large | Tall |
|--------------------------|-----------------|-----------------|-----------------|-----------------|-----------------|
| Code | 7004-0401-000 | 7004-0402-000 | 7004-0403-000 | 7004-0404-000 | 7004-0405-000 |
| 3 Cell Code | 7004-0411-000 | 7004-0412-000 | 7004-0413-000 | 7004-0414-000 | 7004-0415-000 |
| Length | 375mm | 450mm | 550mm | 560mm | 560mm |
| Width | 410mm | 490mm | 580mm | 800mm | 1000mm |
| Height | 100mm | 180mm | 240mm | 240mm | 240mm |
| Depth of W | 80mm | 150mm | 200mm | 200mm | 200mm |
| Width within W | 75mm | 90mm | 125mm | 125mm | 125mm |
| Width at top of W | 75mm | 90mm | 100mm | 100mm | 100mm |
| Cover | 7004-0401-000-C | 7004-0402-000-C | 7004-0403-000-C | 7004-0404-000-C | 7004-0405-000-C |



Positioning Overviews



Snooooooze Star (standard), Twin Side Roll, two Curved Wedges (large), W Leg Trough (medium) and a Banana Cushion (standard)



Snooooooze Star (standard), two Tri-fold Positioning Brackets (large), two Curved Wedges (small), T-Roll (large), Banana Cushion



Snooooooze Star (standard), two Tri-Fold Positioning Brackets (large), two Curved Wedges (small), Curved T-Roll (large), Banana Cushion (standard)



Snooooooze Star (standard), Curved wedges two larges and one small, Twin Roll, Banana Cushion (standard)



Example of accommodation of windswept hips with Snooooooze range: Snooooooze Star (standard), two Tri-Fold Positioning Brackets (large), two Curved Wedges (small), Curved T-Roll (large), Curved Wedge (large), Banana Cushion (standard)

