

medifab

Shaping Better Lives



BINGO

Compact &
supportive
design



For children with complex seating needs, the right stroller can be life-changing, offering safe, supportive positioning, the flexibility to adapt as they grow, and the comfort to enjoy daily adventures.



The BINGO range combines clinical functionality with everyday practicality, giving families the freedom to move with confidence, whether at home, in the park, or off the beaten track.

Functional Benefits



Postural Support Options

Customise your BINGO to suit your child's postural needs, with a wide range of head supports, harnesses, lateral trunk supports, and more.



Storage Basket

The BINGO comes standard with a generous underneath storage basket to hold all the necessities!



Adjustable Handle

Adjust to your preferred height for optimal use.



Removable seat unit

Easily remove the seat unit for storage or travelling.



Tray Table

Add a Therapy Tray for support, play and learning.



Easy-Care Covers

Available in three colours: Turquoise Soft, Berry Soft, and Black.

Designed for easy removal and cleaning to keep your stroller fresh and comfortable.

Forward & Rear Facing

Easily remove and switch the seat unit around so your child can see you or see the world!



Grows with Your Child, Effortlessly

The seat unit of the BINGO Evo is adjustable in seat depth and width, back support height, and lower leg length, meaning that the stroller grows with your child.

BINGO Evo Mini to Size 1 Stroller & Cobra Conversion Kit

As children grow, their supportive equipment should adapt with them, helping reduce stress and avoid unnecessary expense. The BINGO Mini mobility base can be easily converted to fit the size 1 seat unit from the BINGO Evo or BINGO Spex stroller systems, allowing your child to continue using a familiar, comfortable stroller. Simply purchase the conversion kit and a new size 1 seat unit for a smart, cost-effective upgrade to meet your child's changing needs.

Product Codes

BINGO Evo Mini to Size 1 Stroller and Cobra Base Conversion Kit:

2880-8080-086

Removable Seat Unit:

BINGO Evo Size 1: 2881-0011-000

BINGO Spex Size 1: 2885-0011-032



45° Tilt-in-Space



Tilt

The unique tilt mechanism allows for 45° of tilt-in-space in both the forward and rear facing position, without compromising stability. This adjustment allows your child to be better supported whilst out in the world, and is easily changed with the lever, located centrally on the back support.



Recline

The BINGO range features an easy-to-use back support recline function, allowing you to tailor the seat depending on your child's needs and the environment they're in. Whether to rest, or sit up and take on the world, this adjustment is easy to do with the simple turn knob on the rear of the back support.

Compact Folded Size

Remove the seat, and the frame is easy to store in a car boot!
The innovative folding mechanism can be released easily on the push bar.



Comfortable Ride

The BINGO's large, smooth-riding wheels provide a comfortable, protective ride on rough terrain and detach easily with a quick-release mechanism. Its adjustable frame suspension absorbs bumps from all directions and can be fine-tuned without tools for customised comfort.

COBRA High-Low Indoor Base

The COBRA indoor base offers 45° tilt-in-space and gas-assisted height adjustment, enabling your child to access their indoor world with ease. Adjustments can be made while they remain comfortably in their seat, making it simple to use at home or in school. The BINGO seat unit easily detaches from the stroller base and clips securely into the COBRA base for seamless transitions between environments.



Product Codes

Size 0 (Mini): 2880-8720-000

Size 1: 2880-8711-000

Size 2 & 3: 2880-8721-000



Twin Stroller Attachment

A mounted connection that links two strollers, allowing support for two children - one seated in front of the other. This unique, easily detachable connection lets you use the strollers together or separately, offering ultimate flexibility.

Please note: This twin stroller connection is only possible between strollers of the same size, or for connection of the mini and Size 1 together, it is not available for the size 3.

Product Code: 2880-8630-000

Please fill out an additional script form to order the 2nd stroller unit with a unique configuration and accessories.

Please note that this can only be fitted at the time of ordering both stroller units. (Medifab cannot retrofit this attachment).

Attendant Brake

Recommended for purchase with Twin Stroller configuration due to the increased combined weight and inertia when both strollers are connected.

Product Code: 2880-8452-000



BINGO Switch Stroller Frame

The stroller for every adventure

Tackle new activities and take on new adventures with the BINGO Switch, our most versatile stroller frame ever!



Swivelling Front Wheel Attachment:

At the heart of this innovative new design lies an interchangeable front end, a unique feature which eliminates the need to purchase multiple bases. Use the standard wheels for everyday journeys or quickly switch in the third wheel for off-road expeditions, delivering a safe, secure experience across all terrain.

Mini/Size 1: 2884-9571-000

Size 2: 2884-9572-000



Features & Benefits



Transportable

Lightweight and compact folding size for ultimate convenience.



Versatile

Available for BINGO Mini, BINGO Evo and BINGO Spex seating units



Quick Change

Quick release mechanisms make adjusting the frame simple.



Interchangeable Front End

Quickly switch between 3- and 4- wheeled configurations depending on the terrain.



BINGO Switch Sibling Seat Carrier

Extend your adventure, together

The Sibling Seat Carrier is a clever front-wheel accessory for the BINGO Switch stroller frame which is designed to carry a second seat. Compatible with phil&teds® products, including the universal car seat adaptor and the voyager™ double kit™, this innovative add-on is a game-changer for families with both a child with additional needs and a typically developing sibling.

Independent Seat Orientation

The BINGO seat and the extra sibling seat can be turned independently to face either the parent or for forward travel. One child could face forwards and the other could be parent-facing for sibling engagement.

Two Size Options

BINGO Switch Add-On for Sibling Seat:
Size 1 (For Mini, Size 1): 2884-4101-000
Size 2 (For Sizes 2 & 3): 2884-4103-000

Exclusive to BINGO
Switch Frames,
including the BINGO
Switch Mini

Compatible with phil&teds® products,
including the voyager™ double kit™ or
many brands of capsules via
their universal car seat adaptor



(Sibling seat
sold separately).

**Quick-Release
Frame Attachment:**
Easy to connect
and disconnect.

Solid PU Castors for
a smooth ride with
low maintenance.

Accessories

BINGO Evo & Mini

(Includes Switch Models)



Chassis Mount Tray

Size 1: 2880-8058-000

Size 2 & 3: 2880-8059-000

(Size 2 & 3 cannot be used
with BINGO Evo Mini).



Spex 5-Point Harness

2006-6551-000

(Recommended for use with BINGO
Evo Mini, and Size 1 BINGO Evo only).



Wonderseat Padded Shoulder Harness

Size 0: 1409-5660-017



Ankle Huggers

Size 1: 2880-6101-000



Occipital Head Support

Merino/Spacer: 2010-4110-099

Wipe-down: 2010-4110-209



Standard Lateral Head Supports

Merino/Spacer:

Small: 2010-4210-099

Large: 2010-4211-099

Wipe-down:

Small: 2010-4210-209

Large: 2010-4211-209



Extended Lateral Head Supports

Merino/Spacer:

Small: 2010-4310-099

Large: 2010-4311-099

Wipe-down:

Small: 2010-4310-209

Large: 2010-4311-209



Oxygen Bottle Holder

2880-8253-000



Swivel Lock

2880-8440-000



IV Pole

2880-8275-000

Extra Bracket: 2880-8275-050



Foot Support Spacer Block

Size 1: 2880-5321-000




Abductor Pommel

Size 1: 2880-1651-346

Size 2: 2880-1652-346

Size 3: 2880-1653-346

(Size 2 & 3 cannot be used with BINGO Evo Mini).

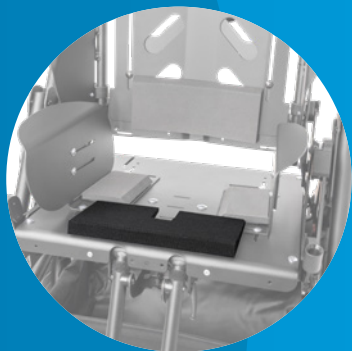

Thigh Wedges

(Recommended for use with the Size 1 BINGO Evo only).

2006-0505-000



Trochanter Wedges

(Recommended for use with the Size 1 BINGO Evo only).

2006-0660-000



Sacral Wedge

(Recommended for use with the Size 1 BINGO Evo only).

2880-2510-060



Wonderseat Hip Belt

4-Point: 1406-5420-027



Hand Strap

2005-6930-000

Accessories

BINGO Evo

(Includes Switch Size 1 & 2)



Standard Laterals

Size 1: 2880-3010-346

Size 2: 2880-3020-346

Size 3: 2880-3030-346



Narrow Straight Lateral Supports

Size 0: 2880-3008-100



Teardrop Lateral Supports

Size 0: 2880-3008-500



Broad Lateral Supports

Size 1: 2880-3008-310



Canopy & Rain Cover

2880-7500-346



Wonderseat Padded Shoulder Harness

Size 1: 1409-5661-017



Wonderseat Hip Belt

2-Point: 1405-6220-027



Foot Support Spacer Block

Size 2 & 3: 2880-5322-000



Waterproof Cover Set - Llamas

Size 1: 2880-9561-121
Size 2: 2880-9562-121
5 Point Harness Cover:
2880-6511-121



Cover Set - Planets

Size 1: 2880-9561-036
Size 2: 2880-9562-036
5 Point Harness Cover:
2880-6511-036



Therapy Tray

Size 1: 2880-7111-000
Size 2 & 3: 2880-7112-000
Infill: 2880-7180-000



Groin Strap

To fit Size 1: 2880-6361-346
To fit Size 2 & 3: 2880-6362-346



Chassis Mount Tray

Size 1: 2880-8058-000
Size 2 & 3: 2880-8059-000



Grab Rail

Size 1: 2880-7251-346
Size 2 & 3: 2880-7252-346



Spacer Fabric Head Support Pads

(For use with Head Support Speakers)

Size 1: 2880-4471-346
Size 2 & 3: 2880-4472-346
Speakers: 2880-8330-000



Centrepont Harness

BINGO 5-point: 2880-6511-346



Chest Shoulder Harness

Size 1: 2880-6751-346
Size 2: 2880-6752-346
Size 3: 2880-6753-346

Accessories

BINGO Evo Mini

(Includes Switch Size 0 / Mini)



Canopy & Rain Cover

2880-7501-346



Grab Rail

Size 1: 2880-7250-346



Standard Laterals

Size 0 (Mini): 2880-3009-346



Broad Laterals

Size 0 (Mini): 2880-3007-310



**Narrow Straight
Laterals**

Size 0 (Mini): 2880-3007-100



Teardrop Laterals

Size 0 (Mini): 2880-3007-500



**Spacer Fabric
Head Support Pads**

Size 0 (Mini): 2880-4470-346



Therapy Tray

Size 1: 2880-7110-000
Infill: 2880-7180-000



**Foot Support
Neoprene Cover**

Size 2 & 3: 2880-5305-000





Here to meet your needs



Needs of parents/caregiver and their child are met through:

- Robust frame construction with multi-directional suspension, providing a more comfortable ride and protection for your child
- Ergonomic, height-adjustable push handle for parents of different heights
- Sun canopy and rain cover option, providing protection for your child when out in the world
- Easy-care upholstery with individual and removable seat, back, head and lateral support sections
- Compact folding for easy storage and transportation, making getting out in the world that little bit easier!

Clinical needs for the child are met through:

- Safe and secure positioning for the user, with comfortable and soft padding
- Adjustable seat unit with 45° tilt in space, back support recline and foot support elevation
- Individually tailored support due to adjustment options for seat depth and width, back support height and lower leg length
- Unrivalled growth capacity with step-less adjustments
- Comprehensive range of scriptable accessories for various postural support and user needs
- Hi-lo Cobra base option for indoor environmental access
- Range of options for transport or support of medical devices (feed pumps, ventilators, supplementary oxygen and suction)

Specifications & Product Order Codes

Size		Code
Bingo Mini	0	2880-0008-000
	1	2881-0008-000
Bingo Evo	2	2882-0008-000
	3	2883-0008-000
Bingo Mini Switch	0	2884-0000-500
Bingo Evo Switch	1	2884-0001-500
	2	2884-0002-500



	Mini (Size 0)	Size 1	Size 2	Size 3
MOBILITY BASE WITH SWIVELLING WHEELS				
Push Handle Height	75–110 cm	75–110 cm	75–110 cm	
Overall Width	63 cm	63 cm	68 cm	
Wheel Diameter Front (Swivel)	19 cm	19 cm	19 cm	
Wheel Diameter Rear	29 cm	29 cm	29 cm	
Tilt in Space (Forward and Rear Facing)	45 °	45 °	45 ° (Forward facing only – Size 3)	
Folded Size Mobility Base (With Rear Wheels on)	80 x 63 x 38 cm	80 x 63 x 38 cm	80 x 68 x 38 cm	
Folded Size Mobility Base (With Rear Wheels Removed)	80 x 63 x 31 cm	80 x 63 x 31 cm	80 x 68 x 31 cm	
Weight Mobility Base (Swivelling Front Wheels)	9.8 kg	8.6 kg	8.7 kg	
Load Capacity Mobility Base	60 kg	60 kg	60 kg	
SEAT UNIT				
Seat Width	13–25 cm	18–30 cm	23–35 cm	23–35 cm
Seat Depth	12–27 cm	16–34 cm	22–40 cm	30–45 cm
Lower Leg Length	16–36 cm	19–40 cm	19–40 cm	19–40 cm
Back Support Height	37–51 cm	45–63 cm	54–75 cm	62–83 cm
Height of Slots for Shoulder Straps	25.5–39.5 cm	28–46 cm	34–55 cm	43–64 cm
Seat to Floor Height (Dependant on Tilt Position of Seat)	50 cm	46 cm	46 cm	46 cm
Backrest Recline	80-180°	90-160°	90-160°	90-135°
Weight of Seat	7.8 kg	6.4 kg	7.4 kg	8.4 kg
Load Capacity of Seat Unit	25 kg	50 kg	50 kg	50 kg
Weight of Seat and Mobility Base Combined	17.6 kg	15 kg	16.1 kg	17.1 kg



medifab

Distributed by Medifab Ltd

solutions@medifab.com

medifab.com

DISTRIBUTED BY:

BINGO Flyer_20250818_AU