

Dreama

Combined pressure-reducing mattress and postural support system, suitable for children and adults. Dreama supports children and adults to achieve supported positioning for lying, helping to reduce muscle tone and contribute to better quality of sleep. Dreama can be used to help reduce the risk of body shape changes and deformity for those who are at risk due to their condition or diagnosis.

Materials

- The cell components that make up the mattress base are made from a pressure-reducing foam, covered in a fire-retardant, water-resistant, 2-way stretch antibacterial fabric. Each cell comes with a comfortable and breathable outer cover as standard, which can be removed easily for washing.
- The majority of the supports available are made from a supportive and comfortable polyurethane foam material. In most cases, these supports are combined with a memory foam pad to provide extra comfort and pressure distribution, in addition to a soft and stretchy towelling outer cover which can be removed and washed.

Cover Options


- All outer covers of the Dreama range are made from a soft and comfortable terry towelling material that assists with breathability and temperature regulation.
- There are no cover options for the Dreama range, however the majority of the supports are waterproof and wipe-clean underneath the outer cover.

Dreama Positioning

- Combined pressure reducing mattress and postural support system, suitable for children and adults.
- Provides supported positioning for at-risk children and adults when lying, helping to reduce muscle tone and contribute to a better quality of sleep.
- Dreama can be used to help reduce the risk of body shape changes and deformity for those who are at risk due to their condition or diagnosis.
- Designed to be very easy to keep clean and hygienic.
- Glide-lock technology means configurations can be changed with minimum disruption to the user. This key feature helps to correctly position posture overnight, whilst providing a secure and stable position for the supports - great for clients who display lots of accessory or uncontrolled movement.



- Contour the flexible aluminium base for a custom lying surface. Suits profiling beds, rolls, wedges.

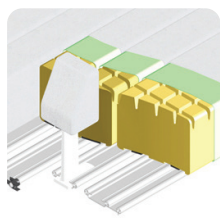


- Quiet with no noise or disturbance when adjusting positions due to Glide-lock supports.



- Quick, simple assembly, and exceptional flexibility with numerous accessories.

Dreama Base & Mattress

- Includes flexible aluminium base, full set of mattress cells, full set of neutral covers and pair of anchor straps.
- Pressure-reducing 15cm deep modular mattress removes the need for overlays.
- Individual covers for each cell means that if soiled, only the covers that need washing can be removed and replaced, saving time and energy.
- Spaces between mattress cushions allow airflow, helping to maintain a comfortable temperature.
- Covers are upholstered in a fire-retardant, water resistant, 2-way stretch, antibacterial fabric for infection control.
- Add locking feet (#7120-6100-000) to mark the position where an accessory should be placed. Locking feet keep accessories in place for users with strong movement.

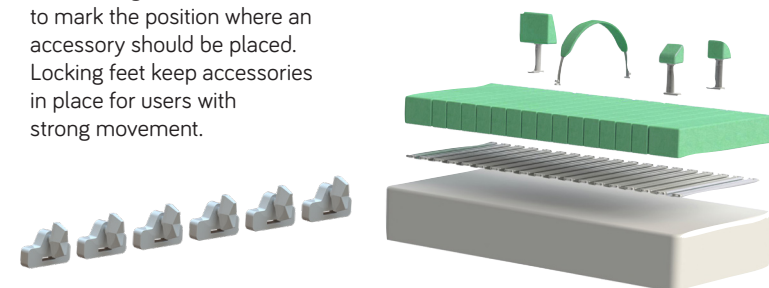


	Single bed
Code	7120-0020-000
Age range	9mths-adult
Overall size	900 x 1900mm
Overall height	170mm
Single mattress cushion	900 x 95 x 150mm
Double mattress cushion	900 x 95 x 150mm
Number of mattress units	15 single, 2 double
Max user weight	110kg

Dreama Positioning Pads

- Made from comfortable and hygienic PU. Positioning Pad Combos include angle-adjustable positioning pads, memory foam pads, and towelling cover.
- Use to position the trunk, arms or legs.
- Slimline Support option is an alternative slimmer support than the standard positioning pads.
- Removable memory foam pad (shown in green) for extra comfort and pressure relief. Ideal for support of areas with bony prominences.
- Additional covers available for all positioning pads in white towelling.


	Pad Combo	Spare Cover	L	W	H
Slimline	7120-3250-500	7120-3280-000	-	108mm	159mm
Small	7120-3251-500	7120-3281-000	102mm	83mm	95mm
Medium	7120-3252-500	7120-3282-000	102mm	83mm	140mm
Large	7120-3253-500	7120-3283-000	203mm	89mm	165mm
X large	7120-3254-500	7120-3284-000	203mm	89mm	267mm

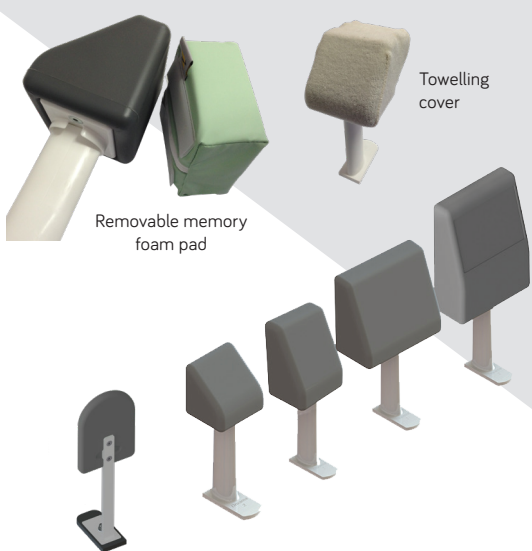

Dreama Abduction Block

- Can be used individually or in pairs to ensure the hips are abducted and the body is straight.
- Ideal for users with increased abductor tone or where positioning into hip abduction is required.


	Small	Medium	Large
Code	7120-1601-000	7120-1602-000	7120-1603-000
Width	130mm	130mm	240mm
Depth	200mm	200mm	300mm
Height	200mm	250mm	300mm
Cover Code	7120-1661-000	7120-1662-000	7120-1663-000

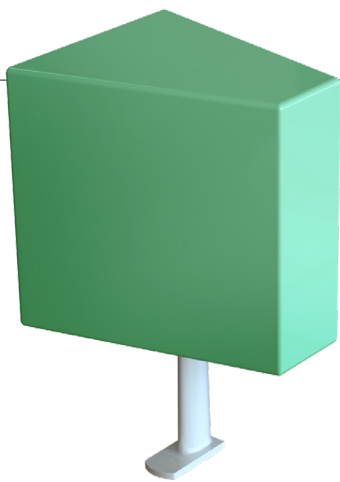

Dreama Knee Blocks

- Comfortably holds the knees in a specified position.

	Small	Large
Code	7120-1751-000	7120-1753-000
Width	77mm	108mm
Depth	102mm	102mm
Height	51mm	102mm
Cover Code	7120-1781-000	7120-1783-000

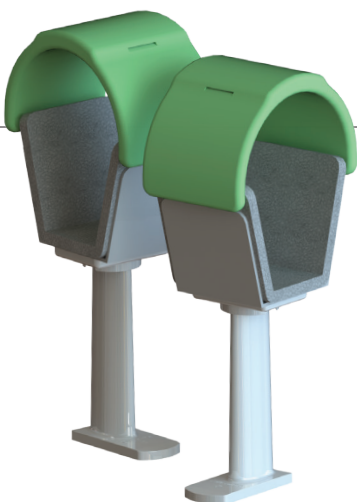

Dreama Side Lying System


- Complete system includes a combination of positioning pads, supportive strap and foam support pad. The straps attach to standard Dreama positioning pads.
- Side Lying Support Strap can be used to secure the upper leg when side-lying. Includes fittings to fit to two standard positioning pads.
- Side Lying Foam Pad is available in different thicknesses to encourage neutral hip positioning of the top hip in side-lying.



	Size 1	Size 2	Size 3
Lying Pad	7120-1851-000	7120-1852-000	7120-1853-000
Pad Cover	7120-1861-000	7120-1862-000	7120-1863-000
Pad Width	345mm	445mm	540mm
Pad Depth	120mm	140mm	159mm
Pad Height	35mm	51mm	76mm
Strap Code	7120-1872-000		
Strap Cover	7120-1882-000		
Strap Dims	140mm wide x 388mm long x 57mm high		

Dreama Positioning Hoop

- A practical solution that helps keep bed covers clear of the body, particularly beneficial for those using medical devices or individuals with sensitive areas.
- Allows movement within a limited range for a comfortable night's sleep.

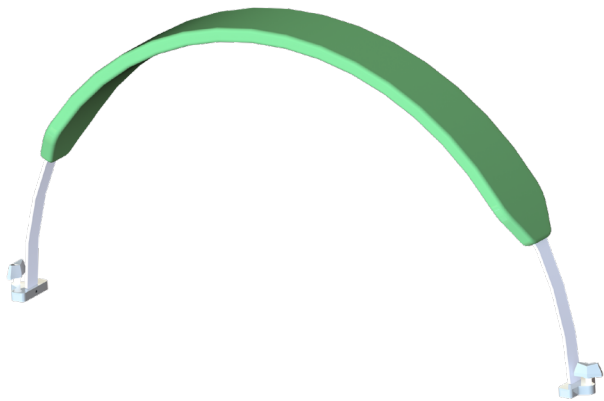


	Small	Medium	Large
Code	7120-1951-000	7120-1952-000	7120-1953-000
Cover	7120-1961-000	7120-1962-000	7120-1963-000
Width	70mm	70mm	70mm
Length	432mm	560mm	700mm
Height	15mm	15mm	15mm

Dreama T-Roll

- The versatile and adjustable Dreama T-Roll was designed to cater to various needs. It provides support for the knees, alleviates abdominal and lumbar pressure, and prevents scissoring and wind-sweeping.
- Each leg can be adjusted independently to accommodate different ranges of motion.
- Removable Medial Knee Support enhances ease of transfer and care management.

Code	7120-1771-000
Width	527mm
Length	108mm
Height	121mm
Cover Code	7120-1761-000

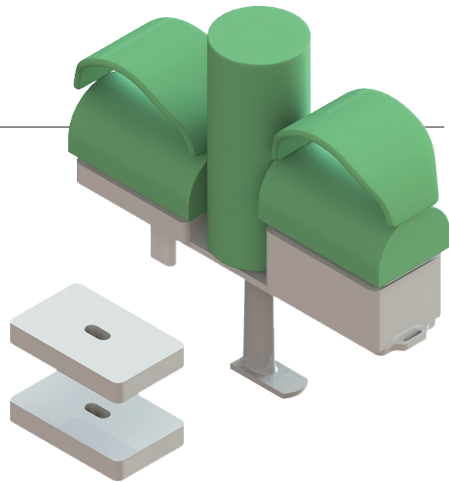

Dreama