

# medifab

Shaping Better Lives



## ENABLE ACC

---

STANDARD & COMPLEX  
LIST EQUIPMENT





# Contents

---

## Simple

|                         |   |
|-------------------------|---|
| Symmetrisleep CoolOver® | 5 |
|-------------------------|---|

## Standard

|                           |   |
|---------------------------|---|
| JAY Ion Cushion           | 7 |
| Spex Vigour Cushion       | 8 |
| Spex Classic Back Support | 9 |

## Complex

|                         |    |
|-------------------------|----|
| QUICKIE 2 Wheelchair    | 11 |
| QUICKIE IRIS Wheelchair | 12 |
| ZIPPIE IRIS Wheelchair  | 13 |

---

## Leaders in Education, Accessible to You

At Medifab, we're more than just providers of high-quality equipment solutions — we're your partners in professional growth and education. Our dedication to continuous learning sets us apart, and we're proud to offer a wide range of educational opportunities designed to support you and your practice. Whether it's through in-person clinics, webinars, or hands-on academies, we're here to help you stay informed, inspired, and equipped to better serve your clients.



### Stay Connected

**Register today** to be the first to know about upcoming education events, webinars and academies.

# Medifab-ACC Simple List Equipment

# Bed and Accessories

CoolOver TR3 is the latest generation of thermal regulation. It provides even superior temperature regulation in that it has:

- Even more thermal capacity – the yarn has been loaded with 50% more of the active ingredient.
- Variable denier weave and hydrophilic properties maximize the wicking; sweat evaporates from the fabric not the body, cooling the fabric which cools the body without chilling.
- Higher wash temperature- 75°C.
- Low temperature tumble dry, high temp for loose sheets.
- No 'whiskering.'

## CoolOver® and thermal comfort

CoolOver® fabrics enhance thermal comfort by working right up close. All-new soft feel Hypo Allergenic fabric contains tiny Omnitherm® microcapsules; which absorb excess heat when the body is hot and release it back when the body cools to maintain a more even temperature for a better night's sleep. Should you turn into the pillow, CoolOver covers have excellent air permeability so reduce the severity of respiratory distress. Inherent free draining characteristics allow any fluids which may otherwise cause redness or rashes to drain away from the skin.

## CoolOver® sheets keep you dry

CoolOver has a 3D 'spacer' knit which allows air to move more freely keeping you cooler as sweat can more easily evaporate and drier as moisture is moved away from the skin.



## Cooler accessories

Not only is CoolOver® available for sheets and covers, it is now used as the standard liner in the Symmetrisleep soft knee blocks and the facing of Symmetrisleep bracket pads as well as the Airmantle and Overmantle. CoolOver covers are available for all Symmetrisleep Lowzone Cushions.

## Implementation and safety checking

It is important to establish a baseline of core temperature (using an aural thermometer) when the individual is well before introducing alternative support surfaces. Remember: Perspiration is not always indicative of a high core temperature.

| PRODUCT                                                         | SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE                  |
|-----------------------------------------------------------------|---------------------|---------------|-----------------------|
| Symmetrisleep CoolOver Single Bed Elasticated Fitted Bed Sheet  | 7015-7120-S00       | ENZA2542      | 106x190cm<br>42x75"   |
| Symmetrisleep CoolOver Single King Elasticated Fitted Bed Sheet | 7015-7120-SK0       | ENZA1430      | 106x210cm<br>42x82.5" |
| Symmetrisleep CoolOver Double Bed Elasticated Fitted Bed Sheet  | 7015-7120-D00       | ENZA695       | 150x215cm<br>59x85"   |
| Symmetrisleep CoolOver Flat Bed Sheet                           | 7015-7120-SFS       | ENZA696       | 210x320cm<br>83x126"  |

# Medifab-ACC Standard List Equipment

# JAY Ion Cushion

The JAY Ion is a multi-layer foam cushion which provides comfort and minimal positioning to users who are uncomfortable in their chair and have a low-pressure risk.

The Ion has a multi-layer foam with a Soft Visco memory foam top layer. For convenience and protection, the JAY Ion features a dual-cover system. The water-resistant inner cover assists with incontinence and moisture management, is easy to clean and includes anti-wicking seam thread and AquaGuard zipper for foam protection.

The machine-washable, stretch outer cover utilises naturally anti-microbial, silver-impregnated X-static® fibre which helps to minimise bacteria build-up (and thus minimise odour). A layer of 3DX spacer fabric promotes airflow for additional skin protection. This cover also helps to reduce shear.



**Max User Weight:** 14"-21" wide: 136kg; 22"-24" wide: 225kg

| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE   | ACC SIC CODE |
|---------------------|---------------|--------|--------------|
| 1705-01-JION1414    | ENZA7847      | 14x14" | EM300        |
| 1705-01-JION1416    | ENZA7848      | 14x16" | EM301        |
| 1705-01-JION1418    | ENZA7849      | 14x18" | EM302        |
| 1705-01-JION1614    | ENZA7850      | 16x14" | EM303        |
| 1705-01-JION1616    | ENZA7851      | 16x16" | EM304        |
| 1705-01-JION1618    | ENZA1514      | 16x18" | EM305        |
| 1705-01-JION1620    | ENZA7852      | 16x20" | EM306        |
| 1705-01-JION1816    | ENZA3148      | 18x16" | EM307        |
| 1705-01-JION1818    | ENZA7853      | 18x18" | EM308        |
| 1705-01-JION1820    | ENZA6609      | 18x20" | EM309        |
| 1705-01-JION2018    | ENZA7854      | 20x18" | EM310        |
| 1705-01-JION2020    | ENZA2149      | 20x20" | EM311        |
| 1705-01-JION2022    | ENZA7855      | 20x22" | EM312        |
| 1705-01-JION2218    | ENZA7856      | 22x18" | EM313        |
| 1705-01-JION2220    | ENZA7857      | 22x20" | EM314        |
| 1705-01-JION2222    | ENZA7858      | 22x22" | EM315        |
| 1705-01-JION2223    | ENZA7859      | 22x23" | EM316        |
| 1705-01-JION2224    | ENZA7860      | 22x24" | EM317        |
| 1705-01-JION2322    | ENZA7861      | 23x22" | EM318        |
| 1705-01-JION2323    | ENZA7862      | 23x23" | EM319        |
| 1705-01-JION2420    | ENZA7863      | 24x20" | EM320        |
| 1705-01-JION2422    | ENZA7864      | 24x22" | EM321        |
| 1705-01-JION2424    | ENZA7865      | 24x24" | EM322        |

# Spex Vigour Cushion

The Vigour Cushion is designed to properly centralise and stabilise the pelvis, hips, and legs to maintain and improve comfort, function, and positioning in all other parts of the body.

Vigour reduces fatigue and pressure with the support and comfort that helps you go the distance. A combination of highly resilient foam and gel-foam that is durable, flexible and comfortable.

**The Vigour Cushion is an excellent choice for individuals who:**

- Require pressure reduction and redistribution.
- Require pelvic stability.
- Need moderate control of posture and upper extremities.
- Need moderate control of abduction or adduction propensities.



**Max User Weight:** 10"-15" wide: 90kg; 16"-20" wide: 160kg

| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE   | ACC SIC CODE |
|---------------------|---------------|--------|--------------|
| 1103-0808-OSC       | ENZA7866      | 8x8"   | EM323        |
| 1103-0810-OSC       | ENZA7867      | 8x10"  | EM324        |
| 1103-1010-OSC       | ENZA7868      | 10x10" | EM325        |
| 1103-1012-OSC       | ENZA7869      | 10x12" | EM326        |
| 1103-1212-OSC       | ENZA7870      | 12x12" | EM327        |
| 1103-1214-OSC       | ENZA7871      | 12x14" | EM328        |
| 1103-1414-OSC       | ENZA6849      | 14x14" | EM329        |
| 1103-1416-OSC       | ENZA7872      | 14x16" | EM330        |
| 1103-1418-OSC       | ENZA3048      | 14x18" | EM331        |
| 1103-1616-OSC       | ENZA4583      | 16x16" | EM332        |
| 1103-1618-OSC       | ENZA1588      | 16x18" | EM333        |
| 1103-1620-OSC       | ENZA7873      | 16x20" | EM334        |
| 1103-1816-OSC       | ENZA7874      | 18x16" | EM335        |
| 1103-1818-OSC       | ENZA379       | 18x18" | EM336        |
| 1103-1820-OSC       | ENZA3098      | 18x20" | EM337        |
| 1103-2020-OSC       | ENZA7875      | 20x20" | EM338        |
| 1103-2022-OSC       | ENZA7876      | 20x22" | EM339        |
| 1103-2220-OSC       | ENZA7877      | 22x20" | EM340        |
| 1103-2222-OSC       | ENZA7878      | 22x22" | EM341        |
| 1103-2418-OSC       | ENZA7879      | 24x18" | EM342        |
| 1103-2420-OSC       | ENZA7880      | 24x20" | EM343        |
| 1103-2422-OSC       | ENZA7881      | 24x22" | EM344        |

# Spex Classic Back Support

Easily solve upper torso needs for most wheelchair users with the Spex Classic. The Spex Classic Back Support incorporates the popular Spex Contouring Cube System to allow on-the-spot customisations for torso comfort and asymmetry needs.

- Effortless customisations
- On-the-spot adaptations for changing postural needs
- Instant results
- Cubes can be removed and added at any time
- Contouring doesn't move or slip out of adjustment
- Additional contouring support wedges can be used to further enhance the degree of contouring



**Standard Hardware:** Quad Mount

**Max User Weight:** Up to 15" wide: 90kg; 16-20" wide: 136kg

| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE   | ACC SIC CODE |
|---------------------|---------------|--------|--------------|
| 1255-1012-OSS       | ENZA7936      | 10x12" | EM426        |
| 1255-1014-OSS       | ENZA5817      | 10x14" | EM427        |
| 1255-1016-OSS       | ENZA7937      | 10x16" | EM428        |
| 1255-1214-OSS       | ENZA5817      | 12x14" | EM429        |
| 1255-1216-OSS       | ENZA7938      | 12x16" | EM430        |
| 1255-1218-OSS       | ENZA7939      | 12x18" | EM431        |
| 1255-1416-OSS       | ENZA5096      | 14x16" | EM432        |
| 1255-1418-OSS       | ENZA5944      | 14x18" | EM433        |
| 1255-1420-OSS       | ENZA7941      | 14x20" | EM434        |
| 1255-1618-OSS       | ENZA4455      | 16x18" | EM435        |
| 1255-1620-OSS       | ENZA385       | 16x20" | EM436        |
| 1255-1623-OSS       | ENZA3351      | 16x23" | EM437        |
| 1255-1818-OSS       | ENZA381       | 18x18" | EM438        |
| 1255-1820-OSS       | ENZA388       | 18x20" | EM439        |
| 1255-1823-OSS       | ENZA390       | 18x23" | EM440        |
| 1255-2018-OSS       | ENZA4487      | 20x18" | EM441        |
| 1255-2020-OSS       | ENZA4142      | 20x20" | EM442        |
| 1255-2023-OSS       | ENZA4404      | 20x23" | EM443        |
| 1258-2220-OSS       | ENZA4645      | 22x20" | EM444        |
| 1258-2420-OSS       | ENZA7942      | 24x20" | EM445        |
| 1258-2620-OSS       | ENZA7696      | 26x20" | EM446        |

# Medifab-ACC Complex List Equipment

# QUICKIE 2 Wheelchair

The modern QUICKIE 2 manual wheelchair series brings advanced aesthetics and technology to the lightweight folding QUICKIE range.

The QUICKIE 2 remains as versatile and modular as before, but is now manufactured out of 7000 Series Aerospace Aluminium which makes the redesigned cross brace more rigid, stronger and lighter.

The stronger folding wheelchair frame results in an increased weight capacity of 136kg, while the overall chair weight is reduced by up to 1kg. The incredibly versatile axle plate along with the modern wheels, caster housings and upholstery provide a great-looking, easy-handling, and even more lightweight folding wheelchair.

Requests for this item **MUST** include:

\* Enable New Zealand's Manual Wheelchair Specification Form & Quote



Accessories ENZ List Code: ENZA3  
Accessories ACC SIC Code: EM153A

Max User Weight: 136kg (standard) and 158kg (heavy duty)

## QUICKIE 2 base unit includes:

- DAB depth and angle-adjustable back support
- 24" lite spoke rear wheels
- Anodised push rims
- Marathon plus tyres
- 5" x 1.5" semi-pneumatic castors
- Single post height-adjustable armrests with full length pads
- 80° swing-away footrest hangers
- Aluminium angle-adjustable footplates
- 16" high 3DX tension adjustable
- Back upholstery
- Integrated push handles
- Anti-tips
- Black Opal frame colour


| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE                 | ACC SIC CODE |
|---------------------|---------------|----------------------|--------------|
| 1571                | ENZA238       | 12-19x12-19"         | EM150        |
| 1571                | ENZA240       | 12-19x20 & 20x12-19" | EM153        |
| 1571-2020-001       | ENZA250       | 20x20"               | EM154        |

# QUICKIE IRIS Wheelchair

The QUICKIE IRIS Manual Tilt Wheelchair is the ultimate tilt-in-space wheelchair, featuring up to 55° of Intelligent Rotation in Space Technology.

It also brings the broadest range of positioning and caregiver options, numerous adjustments, and is one of the lightest tilt-in-space wheelchairs available.

Tilt-in-Space

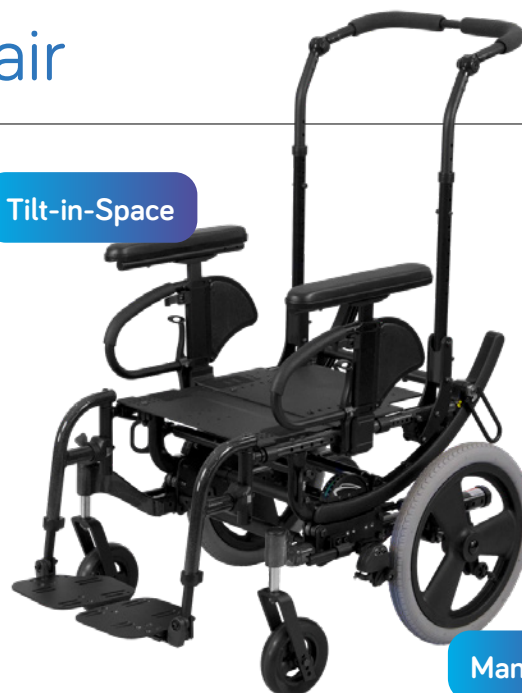

Manual

Requests for this item **MUST** include:

\* Enable New Zealand's Manual Wheelchair Specification Form & Quote



Accessories ENZ List Code: ENZA15  
Accessories ACC SIC Code: EM152A

Max User Weight: 136kg (standard) and 181kg (heavy duty)

## QUICKIE IRIS Medium base unit includes:

- True Fit Growth Program
- Tie downs
- Width adjusting system
- Medium seat frame
- 6" x 1.5" semi-pneumatic front castors
- 16" rear mag wheels; pneumatic
- 80° leg rests hangers with angle-adjustable footplates
- Single post height-adjustable arm supports with full length pads
- Height-adjustable push handle
- Growable seat pan
- Anti-tips
- Black Opal frame colour
- Foot operated tilt and brakes

## QUICKIE IRIS Long base unit includes:

- True Fit Growth Program
- Tie downs
- Width adjusting system
- Long seat frame
- 6" x 1.5" semi-pneumatic front castors
- 16" rear mag wheels; pneumatic
- 80° leg rests hangers with angle-adjustable footplates
- Single post height-adjustable arm supports with full length pads
- Height-adjustable push handle
- Growable seat pan
- Anti-tips
- Black Opal frame colour
- Foot operated tilt and brakes

| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE (width x depth) | ACC SIC CODE |
|---------------------|---------------|----------------------|--------------|
| 1567-00             | ENZA248       | 14-20x16-20"         | EM152        |
| 1567-01             | ENZA259       | 14-20x18-22"         | EM155        |

# ZIPPIE IRIS Wheelchair


The ZIPPIE IRIS uses tilt-in-space technology to rotate the seat frame around the user's centre of gravity.

Available with 40° or 55° rotation range, that can be set to achieve up to 45° or 60° of posterior rotation, the IRIS offers superior positioning capabilities.

This rotational system provides the shortest possible wheel base without compromising stability. It also ensures a smooth and controlled rotation. The ZIPPIE IRIS also comes with the True Fit Program, a complete growth package for free.

Eliminate the need for exposed cables with the Foot Release Tilt Actuator. Never worry again about cable systems getting in the way.

Ability to grow the chair frame as the child's needs change. You can adjust the centre-of-gravity/centre-of-rotation and depth as needed and increase the width with just a few parts.

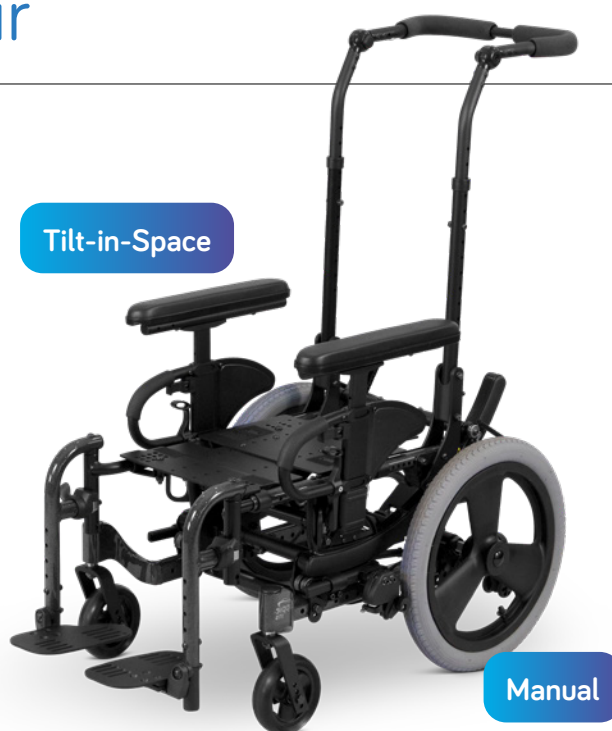

## Requests for this item MUST include:

\* Enable New Zealand's Manual Wheelchair Specification Form & Quote



Accessories ENZ List Code: ENZA14  
Accessories ACC SIC Code: EM151A

Max User Weight: 102kg

## ZIPPIE IRIS base unit includes:

- True Fit Growth Program
- Tie downs
- Width adjusting system
- Short, medium and long seat frame
- 5" x 1.5" semi-pneumatic front castors
- 16" rear mag wheels; pneumatic
- 90° leg rests hangers with angle-adjustable footplates
- Offset hanger mounts
- Single post height-adjustable arm supports with full length pads
- Height adjustable push handle
- Growable seat pan
- Anti-tips
- Black Opal frame colour
- Foot operated tilt
- Foot operated brakes

| SUPPLIER PRODUCT ID | ENZ LIST CODE | SIZE (width x depth) | ACC SIC CODE |
|---------------------|---------------|----------------------|--------------|
| 1562                | ENZA241       | 12-19x12-19"         | EM151        |
| 1562                | ENZA5765      | 10-14x10-18"         | EM151        |

# medifab

Distributed by Medifab Ltd

[solutions@medifab.com](mailto:solutions@medifab.com)

[medifab.com](http://medifab.com)

