

Spex Classic Back Support

Easily solve upper torso needs for most wheelchair users with the Spex exceptional upholstery section. Incorporates the Spex cube contouring cell system to allow on-the-spot customisations for torso comfort and asymmetry needs.

CLINICAL

Designed for the individual who:

- Requires a comfortable back support to reduce torso pain
- Requires spinal alignment and support
- Has skin breakdown/pressure care needs.
- Has bony or protruding skeletal prominences.
- Have changing postural needs
- Pertains a maladaptive or asymmetric posture

FUNCTIONAL

- Comfortable pressure reducing foam surface layer to absorb bony prominences
- Cube contouring cell system provides configurations for a perfect fit and shape for the client
- Customisation to contouring are performed on-the-spot
- Created as the solution-solver for asymmetrical back postures
- Modular componentry means you can 'mix and match' with different lateral pad options and positioning supports
- Offers exceptional back/shoulder contact
- Sustainable choice for issuing to other users
- Dual-layer cover reduces shear forces, and breathable air-mesh on back support sides counters moisture and temperature

RECTANGLE SHAPE

Fits between wheelchair canes



T SHAPE

Sits in front of wheelchair canes






Dual-layer cover reduces **SHEAR FORCES** with air-mesh side panels to counter moisture and temperature

Add laterals as required from our **LARGE MODULAR RANGE**

Strong **VENTILATED SHELL** increases air flow through the back support offering daily comfort

SLOTS & HOLE POSITIONS on shell allow easy installation of head supports, lateral supports and harnesses

CUSHION SIZE CHART

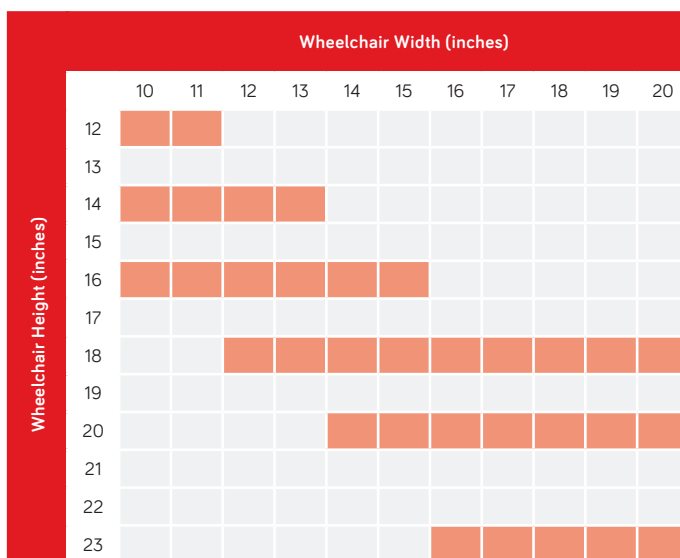

- Red squares are standard sizes.
- Non-standard or custom sizes are available on request.

PRODUCT DETAILS

Width	Capacity	Product Weight
10"-15"	200lbs (90kg)	(16x16)
16"-20"	300lbs (136kg)	1.7kg

- Note: This is the weight of the back support only and does not include the hardware. For hardware weight, refer to page 20.

PRODUCT PERFORMANCE

Postural Control	Red	Red	Red	4
Lateral Stability	Red	Red	Red	5
Postural Asymmetries	Red	Red	Red	5
Changing Needs	Red	Red	Red	5
Pressure Care	Red	Red	Red	5

- Note: This is a guide only and ratings are dependent on correct use of product