

JAY[®]

FUSION[™]

- ≡ All the pressure relief and positioning you expect with added comfort



JAY® FUSION™ CUSHION

≡ The JAY Fusion is an adjustable cushion designed for skin protection and features JAY Flow fluid or air inserts in a package that's comfortable, stable and supportive. It incorporates the stability of the J3 cushion with the pressure management of J2.

Simple to Fit and Order


- Well size based on anthropometric data and sized for cushion width
- Available with air or fluid bladders
- Available in multiple sizes from as small as 14"
- Smart part numbers make ordering easy

Innovative Design Technology


- Multi layered base provides positioning while ensuring comfort
- Soft upper base foam and well fits user comfortably and effectively throughout entire cushion surface
- Incontinent resistant inner cover with Aquaguard™ moisture resistant zipper
- X-static® Anti-Microbial outer cover fends off bacteria and associated smell

Complements Jay Cushion Portfolio

- Base shape and material provides stability and postural support while also providing comfort
- Your choice of Fluid bladder or Air bladder provides required pressure relief
- Jay Fusion provides balance of pressure relief and postural support

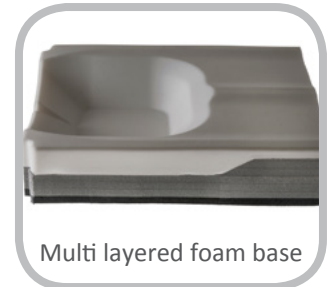

Multi layered foam base

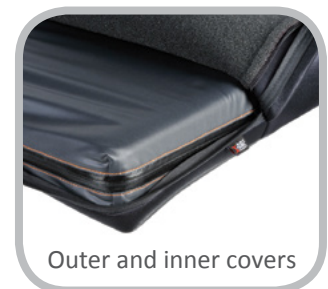

Outer and inner covers



Anti-microbial silver thread

SPECIFICATIONS

JAY FUSION Cushion	
Weight	2.4 kg (16" x 18" with fluid insert) 1.55 kg (16" x 18" with air insert)
Product Width	14" to 24"
Product Depth	14" to 24"
Product Height	4"
Weight Capacity	136 kg (14" to 21" widths) 225 kg (22" to 24" widths)
Outer Cover	X-static®
TGA ARTG #	145061
Warranty	24 months

≡ **JAY**