

JAY[®]

UNION[™]

- ≡ Superior pressure distribution and positioning in a high performance, convenient and anti-microbial cushion




JAY® UNION™

☰ The JAY Union cushion combines multi-layered soft memory foam, a JAY Flow fluid pad and innovative cover technologies to provide outstanding contoured comfort, skin protection and positioning for the client.

Foam and Fluid Composition

- Stable foam base and dual layer Visco™ foam in varying densities create the ideal combination of softness, support, and pressure relief
- JAY Flow™ fluid insert beneath foam fits to each user’s unique pelvic contour
- Contoured base and lateral wedges promote lateral stability
- Pelvic Loading Area provides anterior/posterior stability

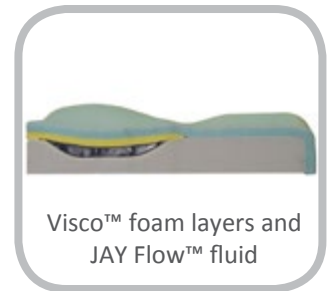

Innovative Cover Technologies


• CONVENIENT

Outer Cover

- Machine washable
- Quick-dry

Inner Cover

- Machine washable
- Aquaguard™ zipper resists moisture
- Anti-wicking seam thread promotes water-resistance to protect foam base

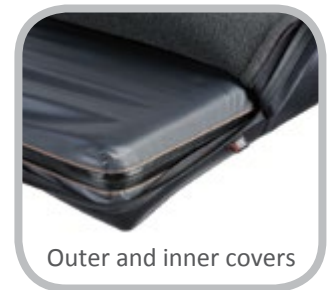

• ANTI-MICROBIAL

Silver-impregnated X-Static® stretch fabric in outer cover

- Naturally inhibits bacteria growth
- Silver thread binds with ammonia and denatured proteins to reduce odor
- Safe non-toxic material



SPECIFICATIONS

JAY Union Back	
Weight	2.4 kg
Product Width	14" – 24"
Product Depth	14" – 24"
Product Height	4.5"
Weight Capacity	14" to 21" widths: 135 kg 22" to 24" widths: 135 kg
Outer Cover	Inner: Water-resistant Outer: Stretch, X-static®
TGA ARTG#	145061
Warranty	24 months

