

Mantaray Back Support

Mantaray offers a solid shell with unique adjustable laterals (or wing tips!) That can be width, height and angle adjusted to give excellent lateral stability. An upholstery wrap offers comfort and secure spine and trunk support.

CLINICAL

Designed for the individual who requires:

- A solid back support to support minimal deformities or stability needs
- Minimal to moderate postural support
- Improved trunk/pelvic positioning whilst maintaining propelling efficiency
- Has mild-moderate torso asymmetries
- Greater functionality as there is minimal hardware between the client and the lateral supports


FUNCTIONAL

- Fully adjustable and repositionable lateral 'wings' means that the therapist/clinician can ensure the right support in the right place for the user
- Slick and non-obtrusive yet provides the right amount of spine and lateral stability for many different postural needs
- Lightweight solid aluminium back support with quick-release hardware
- Wrap upholstery section gives conforming comfort
- Offers exceptional trunk stability
- Minimal contouring, and includes positioning kit for further postural enhancements
- Add head supports and other positioning supports if more support is required
- Approx. Mantaray weight is from 1.5 to 4 kg


CONFORMING COMFORT



COLOURED COVERS AVAILABLE

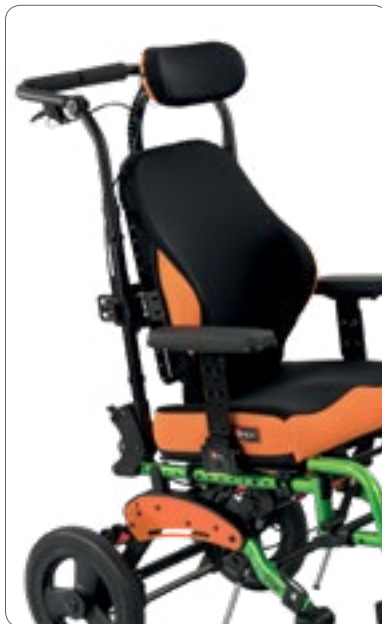


HARDWARE OPTIONS

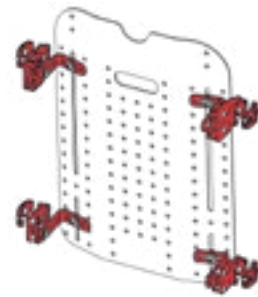

QUAD MOUNT (DEFAULT)

Default hardware with 4 points of attachment for ultimate strength.
#1259-2400-100

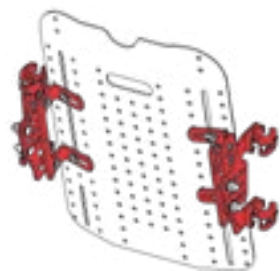

DUAL MOUNT (OPTION)

Increased adjustability, angle and growth.
#1259-2200-100



SINGLE MOUNT (OPTION)

For lower height backs, and as a lightweight option.
#1259-2100-100



Dual-layer cover reduces **SHEAR FORCES** with air-mesh side panels to counter moisture and temperature

Solid back support with integrated **ADJUSTABLE LATERAL SUPPORTS**

Strong **VENTILATED SHELL** increases air flow through the back support offering daily comfort

SLOTS & HOLE POSITIONS on shell allow easy installation of head supports, lateral supports and harnesses

W/C Width	Height	Mantaray Back Support Code	Mantaray Spare Cover Code
12"	14"	1254-1214-000	1254-1214-000-C
"	16"	1254-1216-000	1254-1216-000-C
14"	16"	1254-1416-000	1254-1416-000-C
"	18"	1254-1418-000	1254-1418-000-C
16"	18"	1254-1618-000	1254-1618-000-C
"	20"	1254-1620-000	1254-1620-000-C
17"	18"	1254-1718-000	1254-1718-000-C
"	20"	1254-1720-000	1254-1720-000-C
18"	18"	1254-1818-000	1254-1818-000-C
"	20"	1254-1820-000	1254-1820-000-C
20"	18"	1254-2018-000	1254-2018-000-C
"	20"	1254-2020-000	1254-2020-000-C



MANTARAY WING ADJUSTABILITY

W/C Width	Depth Adjustment	Height Adjustment
Under 16"	50 mm	30 mm
16" and greater	78 mm	40 mm