

Q-SERIES

THE POWER OF QUICKIE



NEW ZEALAND



Advanced engineering. Customised for you.



Q300



Q500/Q500 Mini



Q700



Q700 UP

ADVANCED ENGINEERING. CUSTOMISED FOR YOU.

QUICKIE®, SEDEO®, JAY®, WHITMYER® and SWITCH-IT® bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.



THE POWER OF QUICKIE

A complete power range from basic mobility to complex rehab.
Whatever you need - without compromise.

CLINICAL NEED



Q300

SEATING
SEDEO LITE

BASES
MWD



Q500 / Q500 Mini

SEATING
SEDEO PRO
SEDEO LITE

BASES
MWD



Q700

SEATING
SEDEO ERGO
SEDEO PRO

BASES
MWD
FWD
RWD



Q700-UP

SEATING
SEDEO ERGO UP

BASES
MWD
FWD

□ RWD: Rear-Wheel Drive □ MWD: Mid-Wheel Drive □ FWD: Front-Wheel Drive

BASE PERFORMANCE

SUMMARY

THE NEW QUICKIE Q-SERIES POWER RANGE: A COMPLETE POWER PRODUCT PORTFOLIO UNDER ONE ROOF.

A wide range of bases, drive tractions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.



SERIES

Q700

Q500

Q300



BASES

Front-Wheel Drive

Mid-Wheel Drive

Rear-Wheel Drive



SEATING SYSTEMS

SEDEO ERGO UP

SEDEO ERGO

SEDEO PRO

SEDEO LITE



CONTROLS

VR2

R-Net

SWITCH-IT Controls

Link-It Mounting System

Remote Seating App

SUMMARY

CHOICE WITHOUT COMPROMISE

The new QUICKIE Q-Series offers full flexibility to create the perfect powerchair for each individual user. Highly flexible configurations, interchangeable parts, quick and easy adjustments.



BEST IN CLASS SEATING

Q-Series offers a range of seating systems to cover any clinical need. A choice of six seating systems: SEDEO ERGO UP, SEDEO ERGO, SEDEO PRO, SEDEO LITE, CAPTAIN SEAT and SLING SEAT. The SEDEO range is fully compatible with JAY cushions and backs and WHITMYER headrests to assure the best seating and positioning.

BEST DRIVE BASES CHOICE IN THE MARKET

From cost effective compact bases to high performance outdoor bases, all with the most extensive choice of drive tractions. An impressive range of base levels, and drive tractions that can accommodate all sort of user preferences, lifestyle needs and pricing/reimbursement requirements.

WITH MOST ADVANCED TECHNOLOGIES

Advanced engineering and innovative technologies applied to improve people's lives: SpiderTrac™ 2.0 and Full Track Suspension, Memory Seating, advanced anti-shear backrests, biometrical repositioning, Remote Seating app (pressure management), SWITCH-IT controls and Link-It mounting systems and advanced R-Net Controls.

EASIER TO LEARN

The common design and the modularity create a simplified user experience, that makes it very easy to learn, manage and prescribe.

EASIER TO SERVICE

The modularity and flexibility of the bases, seatings and controls allows an easy adaptation to different users, making it perfect for refurbishments and allowing a great fleet management.

SERIES OVERVIEW

	Q700-UP	Q700	Q500	Q500 Mini	Q300
Outdoor performance	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>
Indoor performance	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>	<div><div></div><div></div><div></div><div></div><div></div></div>
Suspension	<div><div></div><div></div><div></div><div></div><div></div></div> <div>All Wheels</div>	<div><div></div><div></div><div></div><div></div><div></div></div> <div>All Wheels</div>	<div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Advanced</div>	<div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Advanced</div>	<div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Basic</div>
Motors	4 pole	4 pole	4 pole	2 pole	2 pole
Speed	8,10 kph	8,10 kph	10 kph	10 kph	10 kph
Batteries	60 Ah / 72 Ah	60 Ah / 72 Ah	60 Ah / 72 Ah	40 Ah	40 Ah
Max. user weight	120 kg	160 kg	136 kg	136 kg	136 kg
Drive bases	MWD / FWD	MWD / FWD / RWD	MWD / FWD / RWD	MWD	MWD / RWD
Seating	SEDEO ERGO UP	SEDEO ERGO / SEDEO PRO / CAPTAIN SEAT (MWD 136 kg only)	SEDEO PRO / SEDEO LITE / CAPTAIN SEAT (MWD only)	SEDEO PRO / CAPTAIN SEAT (MWD only)	SEDEO LITE / CAPTAIN SEAT (MWD only)
Electronics	R-Net	R-Net	R-Net	R-Net	VR2
Crash tested	Yes	Yes	Yes	Yes	Yes

CREATED INDIVIDUALLY
FOR YOUR CLINICAL AND
LIFESTYLE NEEDS.



SERIES OVERVIEW

QUICKIE Q-SERIES OVERVIEW

Q700 SERIES: THE ULTIMATE EXPERIENCE

Chairs designed to offer a superb outdoor experience while allowing a good indoor use. A choice of three high performance drive bases with market leading suspension technologies, powerful motors and long lasting batteries to provide the ultimate outdoor performance. Combined with a choice of complex rehab seating systems (SEDEO ERGO UP, SEDEO ERGO & SEDEO PRO) to match users with specific or high clinical needs, and the most advanced electronics to access their world. Including also standing versions.

Q500 SERIES: MAKE YOUR OWN ADVENTURE

Robust outdoor powerchairs with a choice of three drive bases with advanced suspension technologies and powerful motors, for a superior outdoor experience. Combined with a choice of rehab seating systems (SEDEO PRO & SEDEO LITE) to provide great support and positioning to users from mild to medium clinical needs, and the most advanced electronics to access their world.

Q500 MINI: ULTRA COMPACT PERFORMANCE SUITABLE FOR COMPLEX REHAB NEEDS

Ultra compact wheelchair with highly adjustable and customisable seating and electronics, a combination to match users with high clinical needs, but tight indoor environments.

Q300 SERIES: ULTRA COMPACT EVERY DAY PERFORMANCE

Entry level chairs with ultra-compact mid or rear wheel drive bases for great indoor agility with the stability and traction for outdoor navigation for community access.





PERFORMANCE BASES

DRIVE BASES OVERVIEW

The new QUICKIE powerchairs range offers the most complete choice of drive bases on the market, available in different series and with many seating systems.

The choice of a drive base affects how the chair performs in different situations and environments. The right choice for each user will depend on personal preferences, experience, clinical, lifestyle and environmental needs.

	MWD	FWD	RWD
Manoeuvrability	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Stability	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Kerb climbing	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Traction	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Drive comfort	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Ease of driving	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
	MWD	FWD	RWD
Q700 UP	✓	✓	—
Q700	✓	✓	✓
Q500	✓	✓	✓
Q500 Mini	✓	—	—
Q300	✓	—	✓



MID-WHEEL DRIVE BASES (MWD)

The mid-wheel drive bases have the best indoor performance. Their smallest turning circle provides great manoeuvrability. They are very intuitive and easy to drive, ideal for users with reduced hand control, driving with switches and users with poor spatial awareness.



FRONT- WHEEL DRIVE BASES (FWD)

The front-wheel drive bases have great outdoor performance with good indoor manoeuvrability. Knee angle can be positioned at 90°. Maintains control on steep declines.



REAR-WHEEL DRIVE BASES (RWD)

The rear-wheel drive bases are stable and have good kerb climbing ability. Intuitive drive that provides good control with non-proportional inputs.



SEATING SYSTEMS

SEATING OVERVIEW

	SEDEO ERGO UP	SEDEO ERGO	SEDEO PRO	SEDEO LITE
User clinical need	High / Standing	High	Medium to High	Medium to Low
Seat width	400 - 500 mm	400 - 560 mm	380 - 560 mm	385 - 510 mm
Seat depth	400 - 540 mm	400 - 560 mm	400 - 560 mm	410 - 560 mm
Backrest height	560 - 610 mm	560 - 700 mm	440 - 620 mm	460 - 510 mm
Standing function	✓	—	—	—
Lay down position	✓	✓	—	—
Programmable seating	✓	✓	—	—
Powered recline	90° - 175° biomechanical	90°-175° biomechanical	90°-130°/90°-150° biomechanical	85°-120°
Manual recline	—	—	90°-120°	90°-120°
Anti-shear back movement	180 mm	180 mm	100 mm (optional)	—
Pressure release reminders	✓	✓	—	—
Tilt angle	22°	50°	50°	30°/50° *
Lift	200 mm	300 mm	300 mm	300 mm
JAY & WHITMYER compatibility	✓	✓	✓	✓
Remote Seating app	✓	✓	—	—
Max. user weight	120 kg	160 kg	160 kg	136 kg
Crash Tested	Yes, ISO	Yes, ISO	Yes, ISO	Yes, ISO

* 50° tilt available on Q500 and Q700 Only



SEDEO ERGO UP

High-end biomechanical seating, with many programmable seating positions and standing function. All the benefits of SEDEO ERGO seating with additional standing function. SEDEO ERGO UP provides a standing angle up to 85° with a highly rigid and stiff standing module.



SEDEO ERGO

High-end biomechanical seating, with many programmable seating positions. SEDEO ERGO covers the most demanding clinical needs with 180 mm anti-shear back movement, biometric repositioning function, infinitely adjustable seat and back, easy entry system, 90° knee angle in seated position and pressure relief reminders.




SEDEO PRO

High-end fully adjustable and configurable seating frame with a wide range of comfort cushions and positioning options. Aimed at meeting moderate to complex clinical needs. With a wide range of armrests, armpads, legrests, seat and back cushions, headrests and positioning accessories, SEDEO PRO offers great flexibility to adapt to the individual user needs.



SEDEO LITE

Multi-adjustable seat frame with the option of either SEDEO LITE backrest cushions, Jay Backrests or sling back upholstery. SEDEO LITE is designed for users with medium to low clinical needs that require adjustable seating system. The backrest, armrests and legrests can be adjusted to fit the user. Includes seat frame with seat rail for additional mounting solutions.

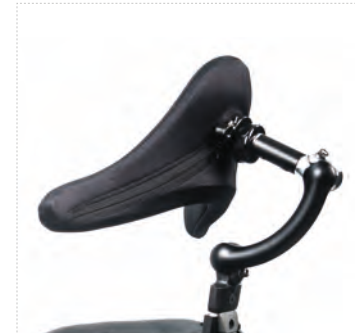

WHITMYER & JAY SEATING

Compatible with all SEDEO seating systems, the JAY family of cushions and backs offers superior clinical seating for pressure management, positioning and comfort. WHITMYER headrests are also compatible with all SEDEO seating systems offering the best head support and positioning.

CONTROLS & CONNECTIVITY

VR2 & R-NET

ADVANCED JOYSTICKS

The QUICKIE Q-Series wheelchairs come with the reliable VR2 control or the superior R-Net advanced joystick. A large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation.

CONNECT

Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick (R-Net only).

BLUETOOTH & IDRIVE

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices (R-Net only).

QUICKIE'S PATENTED ASSIGNABLE BUTTONS

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick (R-Net only).



R-Net
Q500, Q500 Mini,
Q700, Q700 UP




VR2
Q400, Q300,
Q200, Q100

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.

ALTERNATIVE CONTROLS

SWITCH-IT INNOVATIVE CONTROLS

Conveniently access your 5 programmable positions with the Ctrl+5 – a perfect combination with Memory Seating (Memory Seating with SEDEO ERGO and SEDEO ERGO UP only). It allows you to go straight to the function you need at the touch of a button (R-Net only).

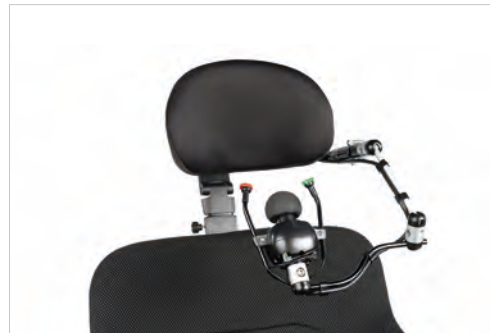

LINK-IT MOUNTING SYSTEM

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It mounting system.







DESIGN AND APPEARANCE

DISTINGUISHED STYLE

With the Q-Series, we've taken style to the next level. With a matt-black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and a contemporary design with a variety of colours to choose from.

Fancy a change?

A choice of contemporary colour inlays to perfectly complement its striking matt-black bodywork. Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



***Q700/Q500:**









***Q500 Mini:**



***Q300:**



* The number of colours available may differ depending on the seat and base combination.

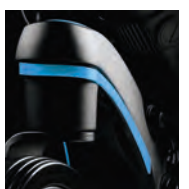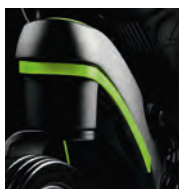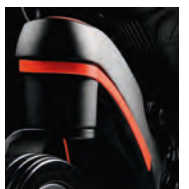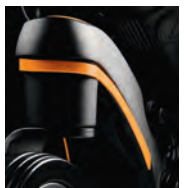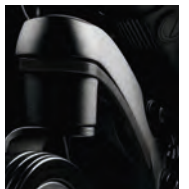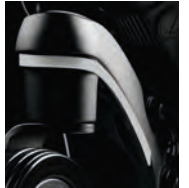



medifab®

Distributed in New Zealand by Medifab Ltd
Head Office: 32 Detroit Drive, Rolleston 7675
solutions@medifab.com
www.medifab.co.nz
0800 543 343

© 2019 Q-Series Range Brochure Au Rev A

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Medifab with any queries you may have. Pictures and illustrations may be overseas models and some seating and accessories not available in New Zealand.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.com.au.

