


Mounting Hardware Options

AXYS

AXYS is a compact mounting system that features a black chrome finish combined with a ½" steel square rod. It was designed for rigorous applications, ease of adjustment, and simple modifications.


Height	Depth	Pitch
16"	6"	30°

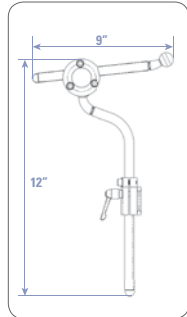

ONYX

ONYX is a compact mounting solution that features a black chrome finish combined with a small profile 7/8" steel rod. It was designed for rigorous applications, ease of adjustment, and simple modifications.


Adult			Paediatric		
Height	Depth	Pitch	Height	Depth	Pitch
11"	Adj 1" to 8"	30°	8"	Adj 1" to 6"	30°

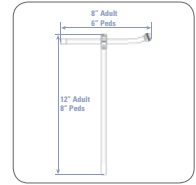

LINX

LINX is a contemporary styled mounting hardware that combines two simple, unique points of adjustment which allows full anterior/posterior fluid movement, creating an infinite range of support surface placement. It features satin black tubing with contrasting polished aluminium clamps and knurled mounting ball ensures stay-in-place adjustments.


Height	Depth	Pitch
7"	8.5" forward, ¼" rearward	30°

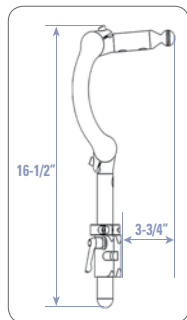

PRO

The PRO Mounting System is a proven hardware option for demanding applications. This mounting system features 1.25" aluminium extrusion and is ideal for use when client strength or strong tonal patterns are an issue.


Adult			Paediatric		
Height	Depth	Pitch	Height	Depth	Pitch
11"	Adj 1.5" forward, 6" rearward	30°	6"	Adj 1.5" forward, 6" rearward	45°

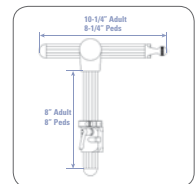

COBRA XTRA

The COBRA XTRA Mounting System incorporates a sleek and unobtrusive design that eliminates protrusions towards the rear of the chair. With superior strength and positioning range, the COBRA XTRA Mounting System is available in Small (Early Intervention and Paediatric), Medium (Paediatric), and Large (Adult). The components are modular and easily changed in the field to enable exact positioning of the head support pad – even if the user's head position changes from the time of the initial prescription.


Adult			Paediatric		
Height	Depth	Pitch	Height	Depth	Pitch
6.5"	6.5" forward, 2" rearward	30°	3.5"	5" forward, 1.5" rearward	30°

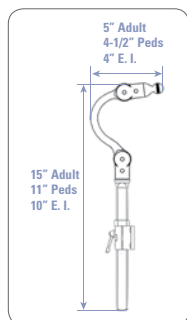

COBRA XTRA FLIP BACK

The COBRA XTRA Flip Back Mounting System is operated by a single push-button on the side of the bracket and neatly flips the assembly back behind the wheelchair back. With superior strength and positioning range, the COBRA XTRA Flip Back Mounting System is available in Small (Early Intervention and Paediatric), Medium (Paediatric), and Large (Adult). The components are modular and easily changed in the field to enable exact positioning of the head support pad – even if the user's head position changes from the time of the initial prescription.


Adult			Paediatric		
Height	Depth	Pitch	Height	Depth	Pitch
6.5"	6.5" forward, 2" rearward	30°	3.5"	5" forward, 1.5" rearward	30°

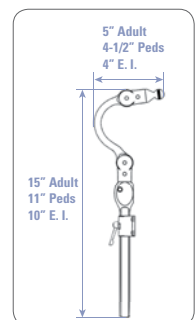

Custom Mounting Options

FLEX INTERFACE BRACKET

The Flex Interface Bracket provides upward flex, travel, and return when the head support is extended against. This dynamic motion assists in maintaining head support contact throughout the user's extensor pattern. The Flex Interface Bracket utilizes elastomeric bushings engineered to supply the correct resistance-to-force while allowing the head support to return to the adjusted position. Used between the head support mount and contact surface, the Flex Interface Bracket can be used with most Whitmyer head support accessories.

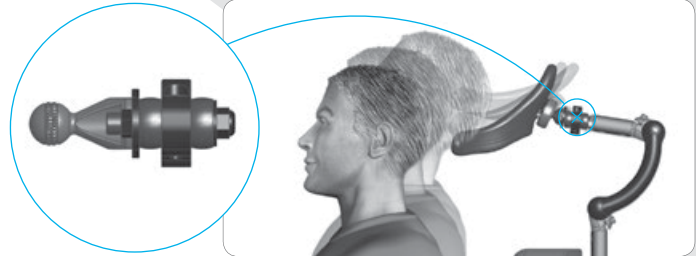

Part Number


M9450

UNIVERSAL SWING-AWAY MOUNT

Whitmyer detach brackets attach to our universal swing-away mount. When attached to the swing-away mount, the assembly becomes a combination detach/swing-away mounting bracket.

Part Number

M50000

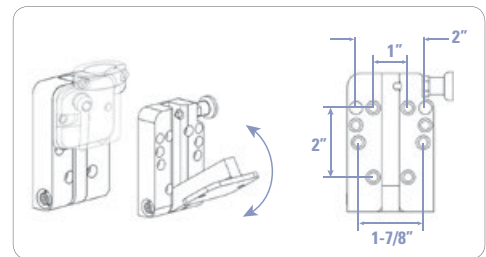


SLIDE MOUNT

This slotted plastic plate allows the head support to be positioned in a more perpendicular fashion behind the head of individuals who present with a significant amount of asymmetry. The plate bolts into the chair's existing hole pattern and allows 4" of movement left or right. Perfect for trying off-set head support locations, adjusting for misplaced T-nut locations and relocating the head support to be more in line with the direction of applied forces.

Part Number

M7010

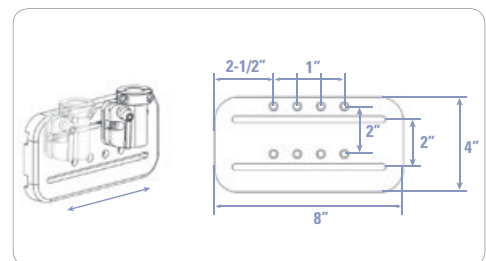


PRO OFF-SET MOUNTING SYSTEM

PRO system off-set bracket provides 2-1/2" of rearward adjustment to the horizontal bar. Flipping the horizontal bar will also allow for 2-1/2" of forward adjustment.

Part Number

PRC7070

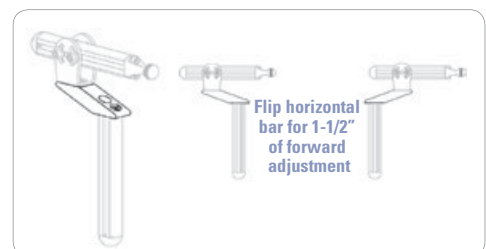

ONYX REARWARD OFF-SET VERTICAL ROD

ONYX off-set vertical rod provides 2-1/2" of rearward/forward adjustment to the horizontal rod.

Part Number

OX6030

