



spex

Shaping Better Lives

POWER
RECLINE
ADAPTER
MOUNT

The Spex Power Recline Adapter Mount provides an off-the-shelf hardware solution for mounting Spex back supports to selected models of power wheelchair systems.

A close-up photograph of a silver metal adapter mount. The mount is rectangular with rounded corners and features the 'spex' logo in the center. It has several mounting holes and slots, including two vertical slots with 'ro' markings. The mount is attached to a black metal frame of a wheelchair. In the bottom left corner, the word 'QUANTUM' is visible on a black component.

spex

The Power Recline Adapter Mount provides an off-the-shelf hardware solution for mounting Spex back supports to selected models of power wheelchair systems.

This versatile Adapter Mount plate empowers clinicians and power wheelchair users by providing an opportunity to select from the Spex seating range to suit individual needs.

Featuring a slim profile and built-in height adjustment, this Adapter Mount maintains the functionality of the back support, while providing a compatible solution for attaching to power wheelchairs.


Compatibility has been tailored to suit two commonly used power wheelchair systems:

- Quantum TRU-Balance 3
- Sedeo Pro (using the J3 Fixed Back Adapter Pro)


Spex has also developed a mounting system for use with the Levo C³ Powered Standing Wheelchair. Regular back support mounting options are also available where power wheelchairs are fitted with back support canes.



Compatible with
Quantum
TRU-Balance 3

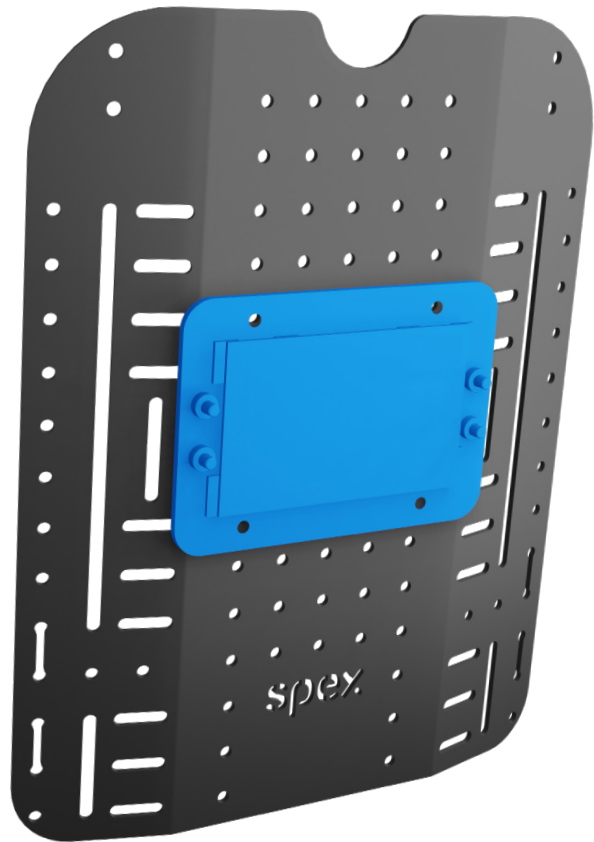

Compatible with
Sedeo Pro
(using the J3 Fixed
Back Adapter Pro)

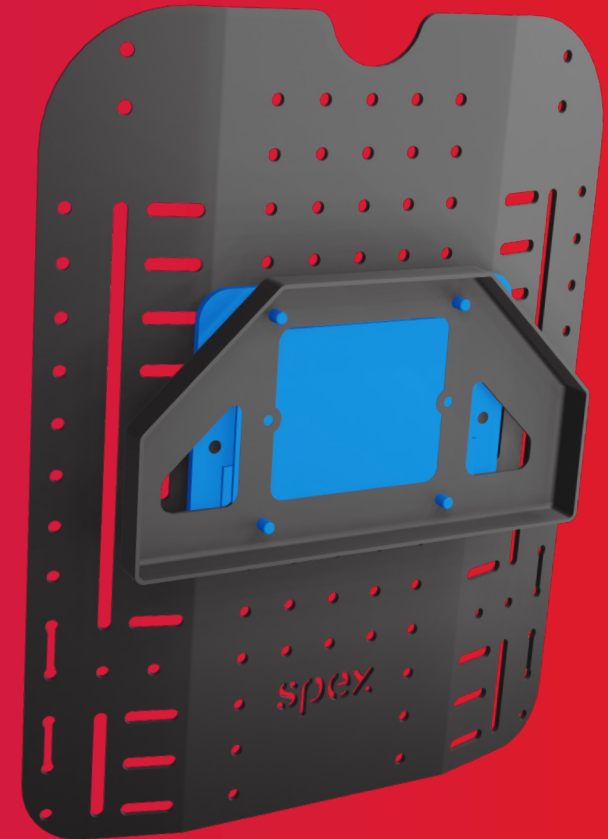

Features & Benefits



Low Depth Profile integrates with Anti-Shear function

The slim depth profile has been designed to work with anti-shear mechanisms used by power wheelchairs for the recline function. This low-depth profile also has the added benefit of minimising seat depth loss.



Positioning Freedom with Height Adjustment

Back supports can be attached to the Adapter Mount at a range of different heights. For fine-tuning adjustments, a further 20mm of continuous height adjustment is available without the need to remove the back support.




Compatibility with Spex Postural Supports


Designed for compatibility with the standard range of Spex postural support accessories, including head supports and lateral supports. Some smaller back support sizes affect compatibility. See Compatibility Guide for further details.

Compatibility with a wide range of Spex Back Supports

The Adapter Mount can be used with a variety of Spex Back Support options, including

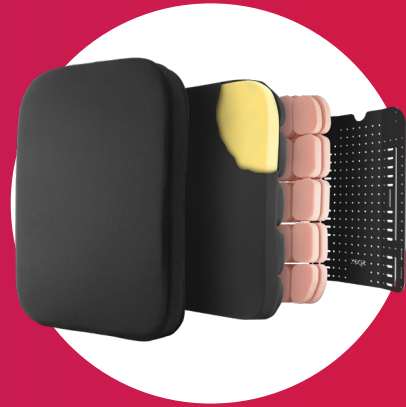

Spex Classic

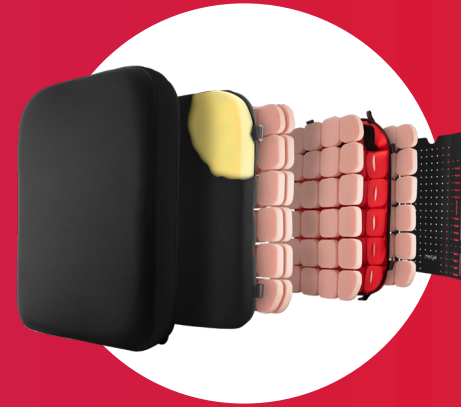

Spex SuperShape



Spex Vigour



Spex Manta



Zygo Mid
Active Contour



Zygo Mid
Deep Contour



XLella Bariatric
Range

Specifications

Power Recline Adapter Mount

Code: 1259-2810-200




Supplied fittings

The side holes are used for mounting to the Quantum TRU-Balance 3 system

XLella Screw Kit for Sedeo Pro

When mounting XLella back supports to the Sedeo Pro system, add this kit to receive longer mounting screws (M6x20 countersink).

Code: 0052-5CSS-SSM6X20BK

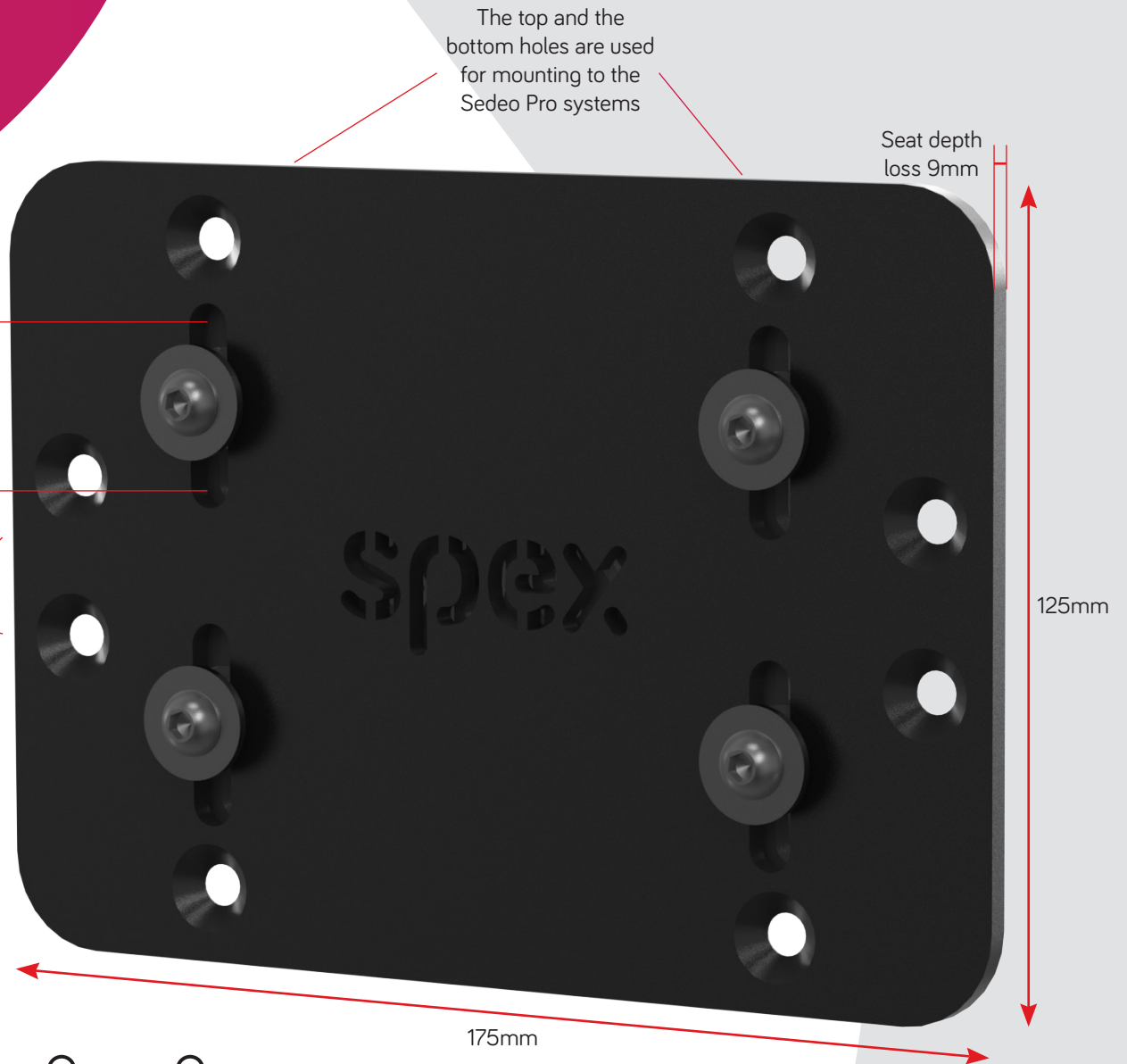



Transit Safety

The Spex Seating System complies with transit safety standards, which evaluate the crash test performance of the seating system independent of a specific model of wheelchair base. The Power Recline Adapter Mount has not been crash tested with a specific wheelchair base. Mounting the Spex Seating System to an appropriate wheelchair base using the Power Recline Adapter Mount, may require a qualified professional to configure and install the seating system.



spex

Manufactured by Spex Ltd

32 Detroit Drive
Rolleston 7675
New Zealand

+64 3 307 9790

solutions@spexseating.com

www.spexseating.com

DISTRIBUTED BY: