

Dreama™ 24hr positioning system by Jenx


At Jenx we have long believed that night-time positioning should be of the same high quality as day time positioning. After all, we probably spend more time in our sleeping position than we do in any other! The importance of 24 hour positioning has become widely recognised over the last few years, as the continuation of good positioning throughout the day and night promotes improved skeletal development and better quality of life. Current thinking is that neglecting proper night-time positioning will take away from the benefits of any postural support given during the waking hours.

An Improved Quality of Life:

The Dreama™ system from Jenx is a key part of any 24hr postural management programme. Recognising that abnormal and posturally destructive tone is present 24/7, the Dreama™ system follows on from support offered by conventional positioning. Utilising a variety of support pads, the system can be configured to accommodate even the most complex of user needs.

Supporting at risk limbs and reducing tonal spasms, the Dreama™ contributes towards more recuperative periods of rest, which in turn has a positive contribution towards minimising the risk of muscular skeletal deformity. Dreama™ can be used as a preventative system or to reduce the risk of further postural deterioration, and is therefore appropriate for any age.

Users who may benefit from Dreama™ include those with the following medical conditions: Cerebral Palsy, Muscular Dystrophy, Spinal Muscular Atrophy, Arthrogryposis, Rheumatoid Arthritis and associated conditions, Post-trauma, Post-surgery, Spina Bifida and Scoliosis.

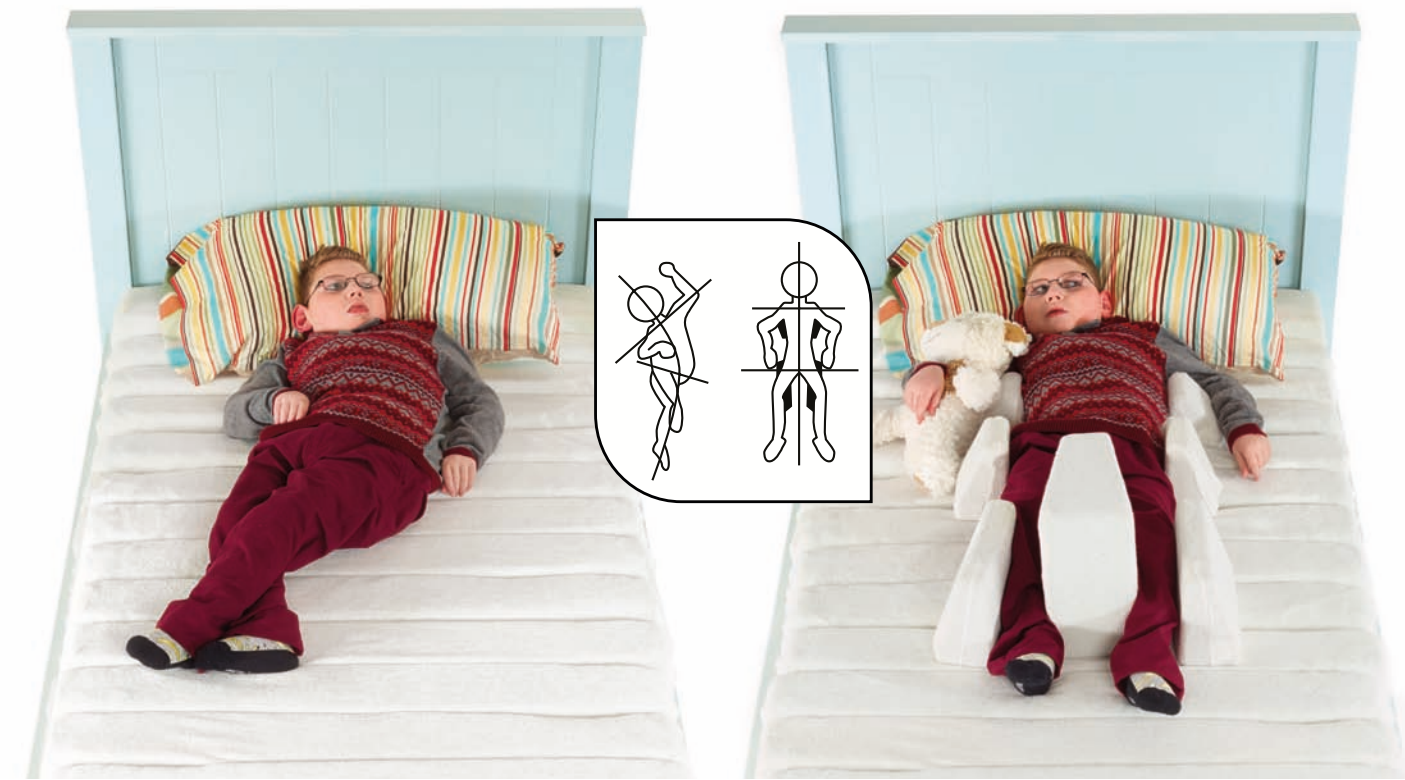

INNOVATIVE FEATURE

Optional Swivel Heads

All pad accessories have the option of a swivel head mechanism. This can be used to change the angle of the support. You can change the angle before, during or after positioning the user.

INNOVATIVE FEATURE

Modular Cushions

Each modular mattress cushion is 15cms (6") deep and upholstered in fire-retardant, water resistant, 2-way stretch, antibacterial fabric for infection control, as well as being breathable. The spaces between cushions allow airflow, helping to maintain a comfortable temperature. The individual cushions enable specific areas to be cleaned easily during the night with minimum disturbance.

Always Supportive

Dreama™ 24hr positioning system - providing valuable rest for everyone

"The Dreama™ sleep system gives Oscar a really good sleeping position overnight, which means he gets a better night's rest and it removes some of our anxiety around ensuring that Oscar maintains the best possible position, but remains comfortable, which means we also get a better night's rest!"

- Julie, Oscar's Mum.

Benefits of Dreama™

Better quality sleep or rest due to comfortable and supportive full depth, pressure reducing mattress.

Easily contours thanks to strong, flexible aluminium base

Allows good ventilation between the body and the mattress, which reduces heat-triggered spasm

Quiet with no noise or disturbance when adjusting positions because of Glide-lock™ supports

Exceptional flexibility for 24 hour support with numerous individual accessories.

Very easy to keep clean and keep hygienic.

Dreama™ by Jenx, is a combined pressure reducing mattress and postural support system.



Showing Dreama™ in side lying

INNOVATIVE FEATURE

Flexible Aluminium Base

The base is fixed to the bed using safety straps and the postural support cushions are then clipped into place.

INNOVATIVE FEATURE

Individual Covers

Individual covers mean any that are soiled can be removed and replaced without the need to remove the user from their bed, minimising sleep disruption and the need for assisted transfers.

INNOVATIVE FEATURE

Glide-lock™ System

The unique Glide-lock™ system makes positioning simple! The base of the pad is inserted in to the channel on the base and the pad is gently pushed in to place. Once in place, simply let go and pushing on the pad won't cause it to move. To remove the pad, simply hold it and carefully slide it out. There are also 'Stops', which can be screwed into place within the channels (they are removable and reusable) either to mark the place a pad should go (making it easier for carers to position a user) or to put behind the pad to prevent users with strong repetitive movements from gradually 'walking' the pad out of place during the night.

INNOVATIVE FEATURE

Customisable and Configurable Positioning

Hundreds of individual specifications can be configured, exactly to the user's needs. The extensive range of positioning options is possible because of the innovative modular mattress. The supports slide easily into channels between the individual mattress cushions. There are 16 available channels for the supports meaning that you can achieve exactly the position required.

INNOVATIVE FEATURE

Custom Profiles

The flexible aluminium base contours to create the lying surface of choice for each user, whether it's achieved by placing the base of Dreama™ on a profiling bed or using rolls, wedges or cushions. This means it can be used for a range of different user requirements in a variety of environments (home, care homes, support centres etc.)

Always Thoughtful Outstanding features and benefits

Pressure Reducing

The main element of the system is an independently-tested pressure-reducing 15cm deep modular mattress. There is no need for any overlays, as using any system on top of a pressure reducing mattress stops that mattress from providing pressure relief.

Comfort

Dreama™ is designed to be comfortable. The towelling covers are soft and stretchy. The pressure reducing mattress gives a supportive base. The support pads mean that the user sleeps in a supported position and can't roll over into an uncomfortable one. Finally, the modular mattress system allows air to circulate, keeping the user cool.

Versatility

Dreama™ offers a full range of positioning opportunities: prone, supine and side lying can all be achieved easily, quickly and quietly.

The individual covers for each pad and mattress cushion mean that you can always see how each accessory is positioned and micro adjust for the perfect amount of individual support.


Quick & Simple set-up

Dreama™ can be assembled within a few minutes, with the mattress cushions clipping into place on the roll out base via integral clips. Configurations can be changed easily without disturbing the user. Supports can be adjusted or the position changed completely with the user still in bed, even asleep! This means less sleep disturbance for everyone!

The flexible base design means Dreama™ can be used on a profiling bed too.

Lifting and Handling

Dreama™ can be set up with the user in bed. This means transfers can occur with some, all, or no supports in place, reducing the need for lifting. It also makes lifting easier as the position can be adapted once the user is on the bed.

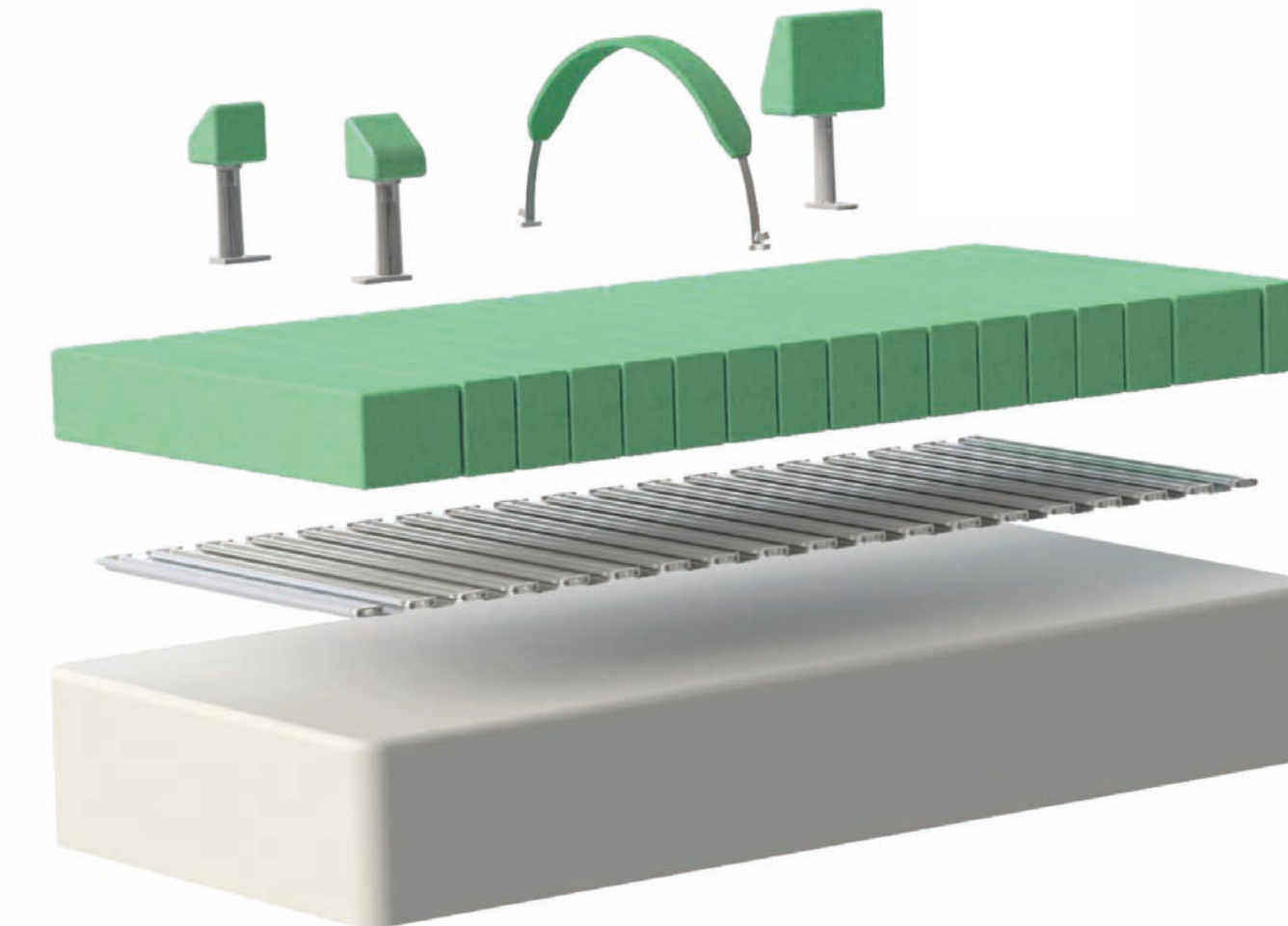


Visit jenx.com for additional product details.

Always Configurable Accessories

The wide range of accessories and optional extras mean the Jenx Dreama™ can be highly customised for every individual's needs. Visit jenx.com for additional accessory details.

Dreama™ Positioning Pads

Made from comfortable and hygienic PU and available in 4 sizes. With or without memory foam pad for extra comfort.

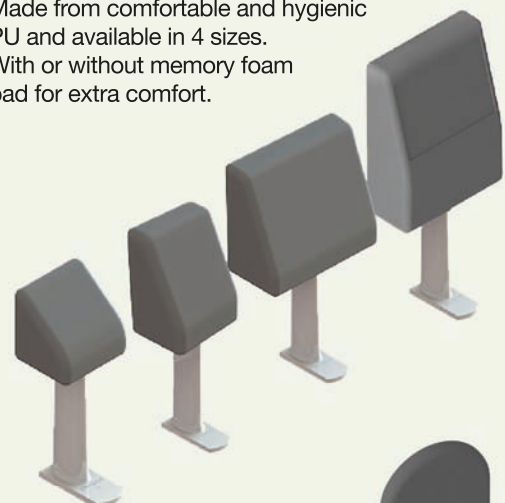

Dreama™ Slim Support

Slimline alternative to the positioning pads. Made from comfortable and hygienic PU. Also available with optional padded cover.



Dreama™ Abduction Block

This can be used individually or in pairs to ensure hips are abducted and the body is straight. Available in 3 sizes.




Dreama™ Knee Cups

Use these to comfortably hold knees in a specified position. Comes in 2 sizes.



Dreama™ T-Roll

This is height adjustable. This helps to support knees, reduce abdominal pressure and prevent scissoring. The sides are height adjustable and the centralommel is removable - a very versatile item.

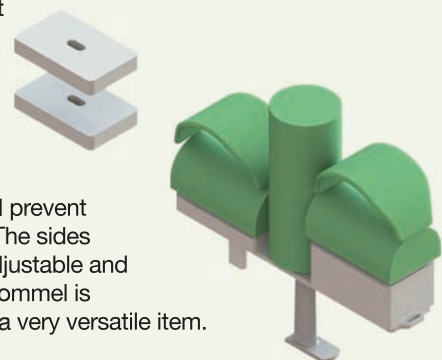

Dreama™ Angle Adjustable Footrest

Ingeniously designed to support feet in supine and prevent dropping.




Dreama™ Positioning Hoops

A brilliant, but functional item, to keep bed clothes away from the body or allow movement within a limited range for a comfortable night's sleep. Available in 4 sizes.




Dreama™ Side Lying System

This complete system includes combinations of positioning pads, straps and Side Lying (abduction) Foam. The straps attach to standard positioning pads with three strap lengths available. The Side Lying Foam is also available in 4 sizes.

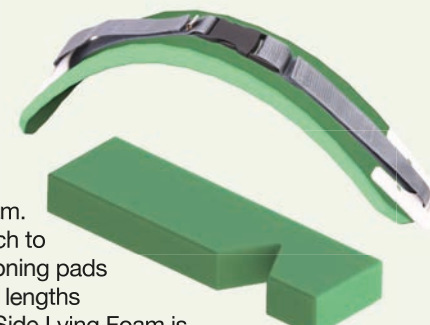


Dreama™ Locking Feet

These useful items are used to mark the position where an accessory should be placed. It keeps accessories in place for users with high movement.

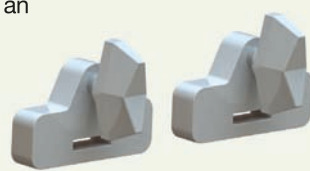

Dreama™ Covers for Support Pads

There are covers available for all Dreama™ accessories, plus various fitted sheet options too. These are provided in white towelling, but can easily be dyed at home to any colour you prefer. See jenx.com for details of all available Dreama™ cover options.

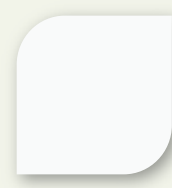

Always Adaptable

Hundreds of individual specifications can be configured with the wide range of accessories



Jenx Dreama™ Technical Specifications

	Units	Cot DC06	Single Bed DC03
Age Range (approx) (Dreama™ Sleep System only)	mths/yrs	9mths-4yrs	9mths-Adult
Overall Dreama™ Footprint	mm	560 x 1200	900 x 1900
Height Overall (Inc.Base)	mm	170	170
Single Mattress Cushion Size	mm	560 x 95 x 150	1200 x 95 x 150
Double Mattress Cushion Size	mm	560 x 195 x 150	1200 x 95 x 150
No. of Mattress Units Supplied (single:double)	-	8:2	15:2

Always here

To arrange a demonstration or assessment with the Jenx Dreama™ 24hour positioning system, just get in touch!

My name is:

Please ring me on:

Or email me at:

Jenx Limited, Wardsend Road, Sheffield S6 1RQ, United Kingdom
tel: +44 (0) 114 285 3376
email: export@jenx.com
web: www.jenx.com



Jenx Limited, Wardsend Road, Sheffield S6 1RQ, United Kingdom
tel: +44 (0) 114 285 3376
email: export@jenx.com
www.jenx.com



Dreama™ by Jenx

The Dreama™ by Jenx is a pressure reducing mattress and postural support system combined, suitable for children and adults.



Always Adaptable...

Visit jenx.com for more information.