

KLAXON®



TWIST

Join the Revolution!

Give your day a Twist



08:00

Good morning!
Start your day with Twist!



08:30

Today is training day!
Marta reaches the gym, with Twist hooked in front to move around safely when outside. Hands free.

08:35

Check-up and set-up
With the Klaxon® App you can set the speed and the braking parameters. Check your battery consumption and much more!



09:35

Fencing training
Marta is a fencing athlete. Twist is perfect for her because she can connect and disconnect it in a moment between the matches.



20:30

Flight to London
Twist, stored in its own bag, is extremely practical to carry and can travel by plane with Pietro. Twist is lightweight and ideal to be taken on all your trips!

11:30

Coffee break
Twist engaged on the front to easily overcome the challenges of the urban routes. Pietro can adjust his speed thanks to the intuitive and ergonomic controller.



19:30

Last-minute shopping
Twist is the ideal solution: it is compact and manoeuvrable even in the tightest of indoor environments such as an airport duty-free shop.

18:30

Emotions are joy
Enjoy unforgettable moments as you hold hands with your loved ones at sunset! With Twist you can reach the most romantic places!

16:00

Playing together
Manuel is an unbeatable dad! He spends his free time in exciting speed races with Erika. Twist is agile, practical and fun!



15:30

In search of a catchy glimpse
Massimo is a painter who draws inspiration from the moving colours. With the handlebar Twist goes faster!



A “thousand” functions Twist!

Several possibilities of use

Twist engaged on the front for a relaxed and safe outdoor ride. Choose your speed with the Controller and move with your hands free. Connect the handlebar and go fast in total comfort! The rear-mounted Twist is ideal for indoor use.



Lightweight and manoeuvrable

Twist is agile. It rotates 360° and is ideal for narrow spaces. The integrated ergonomic handle makes it easy to hold and transport. We take care of every detail.

Quick and easy to connect

Twist can be connected and disconnected in seconds and in total autonomy. Twist is the solution for your independent mobility. Discover a new feeling of freedom without giving up the comfort of moving with your wheelchair.



Compact and transportable

Explore new destinations with Twist. Discover the pleasure of travelling like never before! Twist is compact and lightweight and can be carried and transported by car, train or plane.

Stylish, technological and wireless

Connect your driving system via Bluetooth. Control and select the pre-settings directly from the Klaxon® App. Nothing stands between you and your freedom! Twist is cool and elegant, perfect for moving around in style.



Suitable for all types of terrain

Twist overcomes everyday obstacles with ease. Manages climbs efficiently. When installed in the front it lifts the wheelchair's front wheels and allows to move on all type of terrain.

Safe and secure

Twist is made of high-quality materials and meets the requirements of innovation and technology. It is designed to accompany you at all times in complete safety and durability. Integrated LED lights to live the night!



Linking System Klaxon®

The innovative linking system can be mounted on any type of wheelchair, rigid or foldable. Expand the functionality without sacrificing manoeuvrability and lightness.

Technologies and components

Travel bag
Practical and indispensable bag for your travels.

Integrated LED lights
Ensure visibility and safety.


Foldable handlebar
A handlebar that is practical and easy to carry.



Controller
A touch of the lever is all it takes to accelerate and to choose a pre-set speed. Just press the button to brake.

Removable battery
Twist is even lighter and easier to handle thanks to the removable battery.

Ergonomic handle
Facilitates the engagement and its loading into the various means of transport.

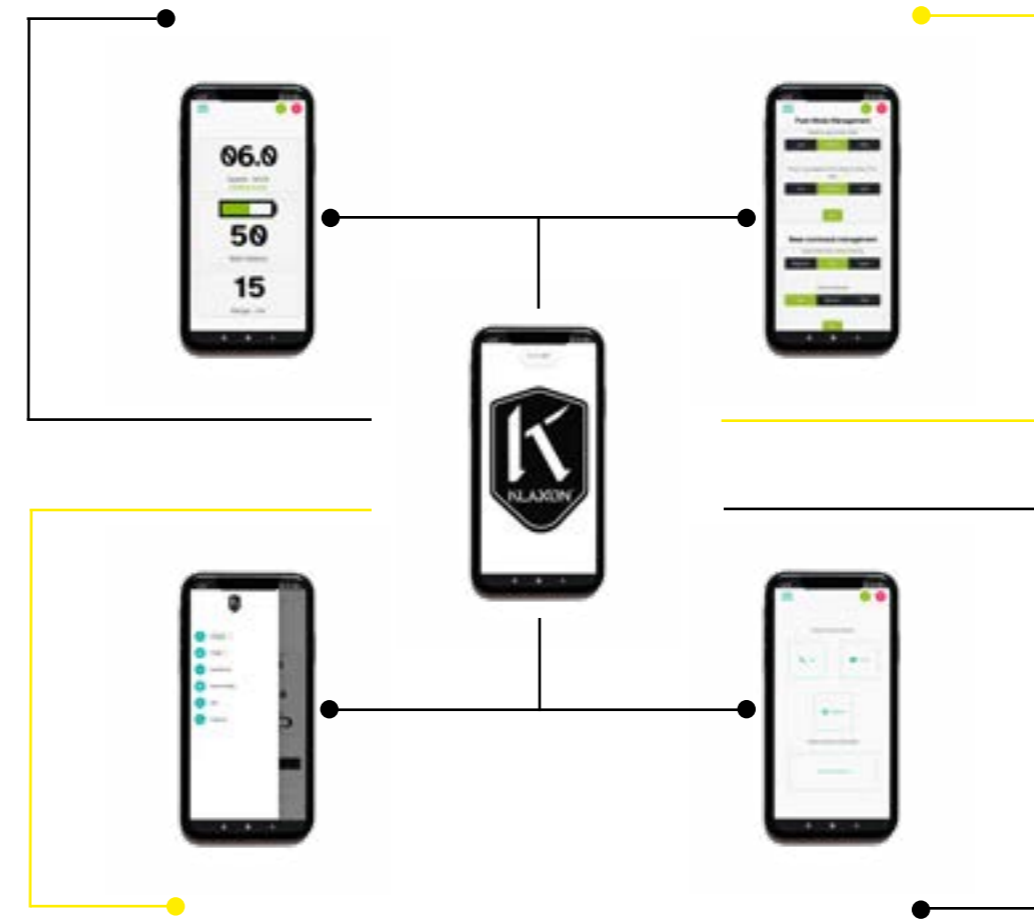

Klaxon® App
Select your favourite driving pre-sets through your smartphone. Stay connected!

Data

Material	Aluminium
Motor	36V x 750W brushless
Wheel	10" smooth tread tyre
Brakes	with Controller: Electronic Brake System with Handlebar: 150 mm Mechanical Brake
Driving speed	with Controller: up to 6 km/h with Handlebar: up to 15 km/h
Battery information	Lithium 36V x 4Ah (weight 1.5 kg) Charging time: approx. 2h
Range	approx. 20 km (at moderate speed, on flat terrain and with user weighing 70 kg)
Weight	7.5 kg (without battery)
Connection	Available options: Central Linking System Klaxon® Lateral Linking System Klaxon® Rear Linking System Klaxon®
Dimensions	Twist: W 15 cm x H 34 cm x D 31 cm with Handlebar: W 45 cm x H 94 cm x D 42 cm
Technologies	Downhill Speed Control (in Controller Mode), Intelligent Cruise Control (in Controller Mode), Electronic Brake System
Included accessories	Battery, Battery Charger, Frontal Light, Rear Light
Included services	Klaxon® App
Optionals	Controller, Handlebar, Additional battery

medifab

Medifab Ltd
22 Detroit Drive
Rolleston 7675
New Zealand

0800 543 343
solutions@medifab.com
www.medifab.com



@medifablives



@medifablives



Medifab



Medifab