

KLAXON®



TWIST join the revolution!

Give your day a Twist

Discover new ways of moving



Safe and secure

Twist is made of high-quality materials and meets the requirements of innovation and technology. It is designed to accompany you at all times in complete safety and durability. Integrated LED lights to live the night!



Suitable for all types of terrain

Twist overcomes everyday obstacles with ease. Manages climbs efficiently. When installed in the front it lifts the wheelchair's front wheels and allows to move on all type of terrain.



Stylish, technological and wireless

Connect your driving system via Bluetooth. Control and select the pre-settings directly from the Klaxon® App. Twist is cool and elegant, perfect for moving around in style.



Quick and easy to connect

Twist can be connected and disconnected in seconds. The innovative Linking System Klaxon® makes it very easy to install Twist on foldable and rigid wheelchair.



Lightweight and manoeuvrable

Twist it is light. The engine's thrust smoothes your manoeuvres. It rotates 360° and is ideal for narrow spaces. We take care of every detail.



Compact and transportable

Twist is compact and can be carried and transported by car, train or plane. The integrated ergonomic handle makes it easy to lift it and transport it. Travelling made easy!

Technologies and components



Controller

A touch of the lever is all it takes to accelerate and to choose a pre-set speed. Just press the button to stop.

Linking System Klaxon®

Thanks to this patented Linking System, Twist can be installed on the front for outdoor and indoor activities. Engaged on the rear it is ideal for indoor mobility.

Integrated LED lights


The two led lights ensure visibility and safety. To live the night!

Ergonomic handles

The handle of the rear linking system is designed to ease the linking/unlinking. The integrated handle facilitates its transportation.

Travel bag

The travel bag allows you to store your Twist in the luggage compartment of your car or in the overhead compartment of an aircraft.

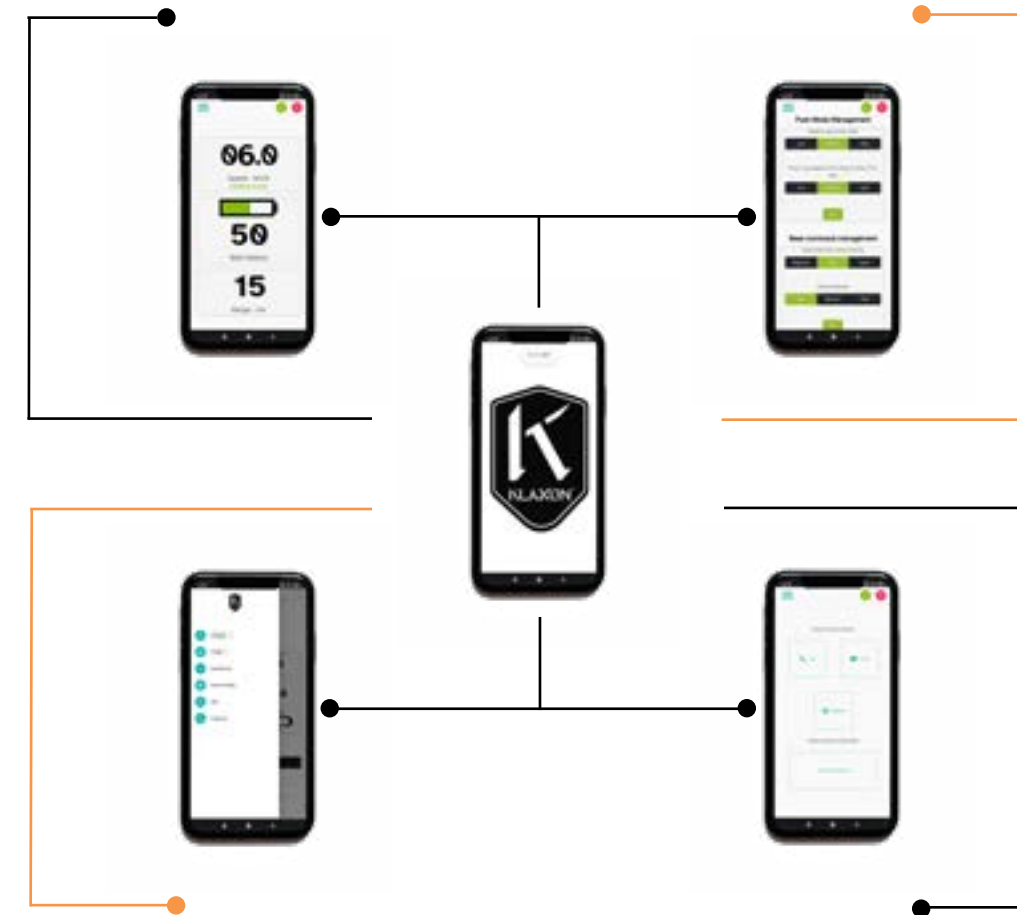

Klaxon® App

Select your favourite driving pre-sets through your smartphone. Stay connected!

Data

Material	Aluminium
Motor	36V x 750W brushless
Wheel	10" smooth tread tyre
Brakes	with Controller: Electronic Brake System
Driving speed	with Controller: up to 6 km/h
Battery information	Lithium 36V x 4Ah (weight 1.5 kg) Charging time: approx. 2h
Range	approx. 20 km (at moderate speed, on flat terrain and with user weighing 70 kg)
Weight	7.4 kg (without battery)
Connection	Available options: Central Linking System Klaxon® Lateral Linking System Klaxon® Rear Linking System Klaxon®
Dimensions	Twist: W 15 cm x H 34 cm x D 31 cm
Technologies	Downhill Speed Control (in Controller Mode), Intelligent Cruise Control (in Controller Mode), Electronic Brake System
Included accessories	Battery, Battery Charger, Frontal Light, Rear Light
Included services	Klaxon® App
Optionals	Additional battery

“Come on, let’s Twist again!”



Let’s Twist Again - Chubby Checker

medifab

Medifab Ltd
22 Detroit Drive,
Rolleston 7675 NZ
0800 543 343

solutions@medifab.com
www.medifab.com

