

QUICKIE® 2



GREAT MADE
Greater!



GREAT MADE Greater!

MODERN TECHNOLOGY LIGHTER WEIGHT

The modern QUICKIE 2 brings advanced aesthetics and technology to the folding ultralight wheelchair. 7000 Series Aerospace Aluminum makes the re-designed crossbrace more **rigid, stronger, and lighter**. The stronger frame results in an increased weight capacity of **136kgs**, while the overall weight is reduced.

The incredibly versatile axle plate along with the modern wheels, caster housing and upholstery, provide you with a great-looking, easy-handling, and even more lightweight folding chair.



Q-FIT™ TECHNOLOGY

The QUICKIE 2 uses Q-Fit Technology to create a tight fit between components for a snug, out-of-the-box feel over the lifetime of the wheelchair. Living springs in the side guards and padded swing-away arms eliminate wiggle and rattling, while reversible seat saddles lock the seat rails in place for a smooth, rigid ride.

ENDLESS OPTIONS CUSTOM FIT

The QUICKIE 2 is customised to your individual angles and measurements. The modular frame and interchangeable components can change with you as your body and environment change. Choose from over **1,000 configurations** of components.

You can also choose the hemi frame, depth-adjustable back, angle-adjustable back, or the QUICKIE Xtender™ Power Assist for extra push.

With more selection of features and more customised measurements, the QUICKIE 2 is built to be an extension of you.



QUICKIE® 2

Made Lighter.....Stronger.....
GREATER!

Quickie 2® is now lighter, accommodates a higher weight and incorporates all the benefits the Quickie brand is known for – strength, durability and style.

The modern Quickie® 2 brings advanced aesthetics and technology to the folding ultra lightweight manual wheelchair.

7000 Series Aerospace Aluminum makes the redesigned crossbrace more rigid, stronger, and lighter. The stronger frame results in an increased weight capacity of 136kgs, while the overall chair weight is reduced. The incredibly versatile axle plate along with the modern wheels, caster housings, and upholstery provide you with a great-looking, easy-handling, and even more lightweight folding wheelchair.




EXO TENSION ADJUSTABLE BACK UPHOLSTERY

Modern, lightweight, and comfortable. Offered with or without tension adjustability, to provide optimal positioning and simple adjustment.



VERSATILE AXLE PLATE

The reversible axle plate allows for exceptional adjustability and options, with fewer parts.

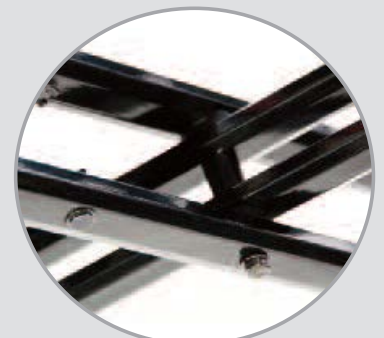

LIGHT & DURABLE CROSSBRACE

Made with 7000 Series Aerospace Aluminum, the diamond-shape crossbrace arms are engineered for strength and greater rigidity.



HEIGHT ADJUSTABLE PUSH HANDLES

Easy height adjustment without added hardware or receivers. The push handle is designed for an easier push, while the folding design allows easy storage.

Q2 LITE & Q2 HP OPTIONS

The QUICKIE 2's minimalist counterpart, the QUICKIE 2 Lite, features an extremely lightweight, modular frame with a rounded rear for unique style and either fixed or swing-away front hangers.

The QUICKIE 2 HP™ or High Performance, was built for an active user looking for the most rigid ride available from a folding frame.



QUICKIE[®] 2

SPECIFICATIONS

Overall Weight
from 11.5kgs

Transport Weight
from 6.2kgs

Weight Capacity
Standard: 136kgs
Heavy Duty: 158kgs

Seat Width
12" to 22"

Seat Depth
12" to 20"

Back Angle Adjustment
Fore (-) 3° to Aft (+) 15°

Front Seat Height
13.5" to 22"

Rear Seat Height
13" to 22"

Hanger Angle
60°, 70°, 80°, 90°

Transit Approved
Yes

Warranty
Limited/Lifetime

medifab

Distributed in New Zealand by Medifab
Phone: 0800 543 343
solutions@medifab.com
www.medifab.co.nz