

Q400[®] SERIES

Advanced Outdoor/Indoor Performance



Q-Series
Powerchairs

Advanced engineering. Personalised for you.



Q400 M
Lite Seat

ADVANCED ENGINEERING. Personalised FOR YOU.

QUICKIE®, SPEX SEATING, SEDEO®, JAY®, WHITMYER® and SWITCH-IT® bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.



THE POWER OF QUICKIE

From basic mobility to complex rehab. A wide range of bases, drive wheel positions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.

CLINICAL NEED



Q300 Mini

SEATING
SEDEO LITE

BASES
MWD



Q400

SEATING
SEDEO LITE

BASES
MWD



Q500 / Q500 Mini

SEATING
SEDEO PRO
SEDEO LITE

BASES
MWD
MWD



Q700

SEATING
SEDEO ERGO
SEDEO PRO

BASES
MWD
RWD
RWD



Q700-UP

SEATING
SEDEO ERGO UP

BASES
MWD
FWD

■ RWD: Rear-Wheel Drive
 ■ MWD: Mid-Wheel Drive
 ■ FWD: Front-Wheel Drive

BASE PERFORMANCE

ADVANCED OUTDOOR AND INDOOR PERFORMANCE

The new QUICKIE Q400[®] offers great outdoor performance with good indoor manoeuvrability and a head turning style. All these combined with Spex Seating or the multi-adjustable SEDEO LITE Seating to cover a wide variety of clinical needs.



PERFORMANCE BASES

Q400 M: Mid-Wheel Drive

SEATING SYSTEMS

SPEX SEATING

SEDEO LITE

JAY & WHITMYER

CONTROLS & CONNECTIVITY

VR2 Control

DISTINGUISHED STYLE

STYLISH DESIGN and choice of colours

PERFORMANCE BASES

Q400 M MID-WHEEL DRIVE BASE

The QUICKIE Q400 M base effortlessly balances great outdoor manoeuvrability with very good indoor performance. Its motors combined with 6 wheels suspension and QUICKIE's anti-pitch technology, and a compact turning radius, will ensure you have a great driving experience.

SUPERB INDOOR MANOEUVRABILITY

Great access through narrow doors, in lifts and around confined spaces, with the compact turning radius and base dimensions. Easily transfer in and out, effortlessly access vehicles and under tables with the lower than ever before seat to floor height from 41.5 cm with powered COG Tilt.



GREAT OUTDOOR PERFORMANCE



ALL-WHEEL SUSPENSION

Enjoy a smooth and stable ride, even over rough terrain, with QUICKIE's patented all-wheel suspension - six independently suspended wheels absorb all those bumps. Climb kerbs of up to 100 mm and have the best traction even when ascending gradients of up to 8 degrees.



ANTI-PITCH TECHNOLOGY

QUICKIE's anti-pitch technology eliminates pitching on slopes, and with all your weight located directly above the drive wheel to maximise traction and drive performance.



ENJOY YOUR TIME OUTSIDE

With reliable 4-pole motors powering speed of 10 kph and 60 Ah batteries to maximise your range, you can always enjoy your time outside with the Q400 M.



SEATING: SEDEO LITE



	SEDEO LITE
User clinical need	Medium
Seat width	385-510mm
Seat depth	410-560mm
Backrest height	460-510mm
Configurable	Medium
Standing function	—
Lay down position	—
Programmable seating	—
Powered recline	85°/120°
Anti-shear back movement	—
Pressure release reminders	—
Tilt angle	30°
Lift	300 mm
SPEX SEATING, JAY & WHITMYER compatibility	✓
Remote Seating app	—
Max. user weight	136 kg



SEDEO LITE: MULTI-ADJUSTABLE SEATING

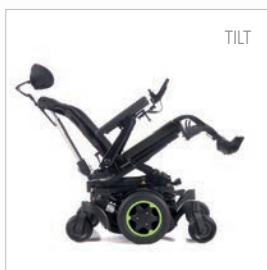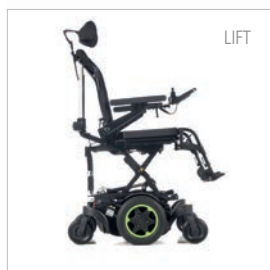

*Image showing available recline option

SLING BACK WITH OPTIONAL COMFORT CUSHIONS

The SEDEO LITE offers height adjustable sling back with adjustable straps, and a range of optional comfort cushions including two seat cushions and four back cushions with adjustable straps and different contours for better comfort and positioning.

JAY & SPEX cushions and backs and are also compatible with SEDEO LITE frame to cover higher clinical needs.

FLIP BACK ARMRESTS

The easy to operate flip back armrest gives you free transfer space. You can choose from a removable or reclining version. All compatible with different armpad options.

LENGTH COMPENSATING SWING AWAY LEGRESTS

The swing away length compensating leg-rest allows you to reach the best position. Also available a wide range of swing away and centre mount leg-rests both in manual and power elevating versions. From simple version to the SEDEO LITE swing away leg-rest with many options.

INTEGRATED HEADREST DESIGN

The integrated headrest mount provides a new standard in adjustment and design. Compatible with SPEX, SEDEO and WHITMYER headrests.

POSITIONING ELEMENTS & INTEGRATED SEAT RAIL

Different positioning elements are available like swing away lateral thoracic supports, thigh supports and knee supports. All of them are easily swing away and detachable. The integrated seat rail offers many options for mounting the different positioning elements.



SPEX, JAY & WHITMYER POSITIONING OPTIONS

A wide range of SPEX cushions, back & head supports, WHITMYER headrests and JAY cushions are available and compatible with the SEDEO LITE frame to meet higher clinical needs and ensure we offer the best seating possible.



WHITMYER HEADRESTS



Created to address a range of user needs from basic support to multiple-point control. WHITMYER headrests come in a variety of pad shapes and sizes for superior head positioning and comfort.

JAY CUSHIONS & BACKS

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly personalised and are available in a large variety to meet your needs.

SPEX CUSHIONS, BACK & HEAD SUPPORTS

Designed to provide exceptional outcomes for all seating requirements. The wide range from Spex includes backs which solve upper torso needs, cushions which provide optimum pelvic foundation and head supports which give adequate head and neck control.

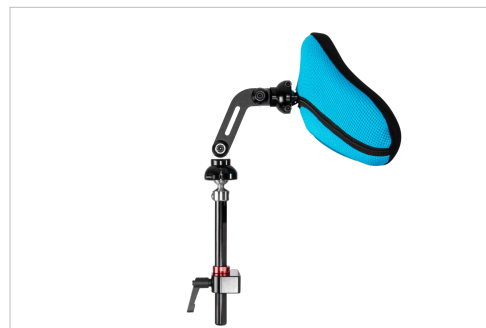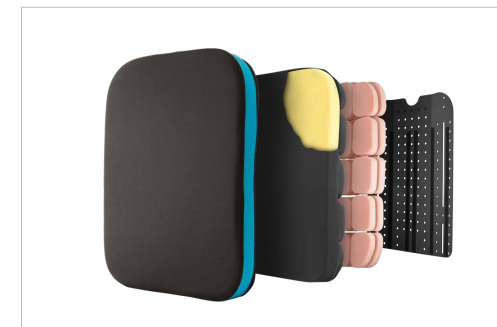

DRIVE CONTROLS

CONTROL AT YOUR FINGERTIPS



SEDEO
LITE

VR2 Joystick

Every Q400 comes with a reliable VR2 control and a simple to control joystick which allows for intuitive operation.



VR2 Joystick

ATTENDANT CONTROL

Allows for a carer to assist with driving or power seat function control.



Attendant Control




DESIGN AND APPEARANCE

DISTINGUISHED STYLE

With the Q400 series, we've taken style to the next level. With a Matte Black chassis finish, hidden wiring and components, it's a sleek and contemporary design.

Fancy a change?

The Q400 series also come with a choice of 4 contemporary colour inlays to perfectly complement its striking Matte Black bodywork. Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.

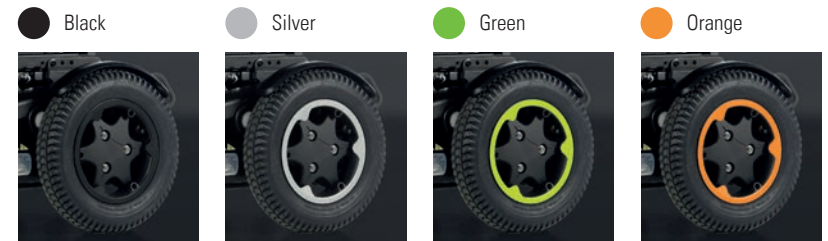

TECHNICAL DATA / Q400 M

	Q400 M w/ SEDEO LITE		Q400 M
		Max. speed	10 kph
Max. user weight	136 kg	Battery size	60 Ah
Seat width	385-510 mm	Electronics	VR2
Seat depth	410-560 mm	Turning radius	560-630 mm
Seat height	415-490 mm	Max reversing width	1100 mm
Backrest height	460-510 mm	Max. safe slope	8°
Total height	Max. 1400 mm	Max. kerb climbing	100 mm
Powered recline	85°-120°	Drive wheels	13"
Tilt angle	30°	Castor wheel options	7"
Antishear back movmt.	–	Colour inlay choices	4 Colour Accents
Seat functions	Tilt, Lift & Tilt	Max. range (ISO 7176-4)	Max. 40 km
		Total width	610-620 mm
		Total length	Max. 1130 mm
		Wheelchair weight	From 130 kg



Sunrise Medical Australia Pty. Ltd.
 11 Daniel Street
 Wetherill Park NSW 2164
 Phone: (02) 9678 6600
 Email: enquiries@sunrisemedical.com.au
www.SunriseMedical.com.au

@ 2020 Q400 Range Brochure AU Rev D

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

Pictures and illustrations may be overseas models and some seating and accessories not available in Australia.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.com.au.

