



# X'PLORE your world!

The Zippie X'CAPE is the first ever wheelchair with XLOCK™ Technology for the performance of a rigid frame and convenience of a folding frame. Designed for active and growing kids, the X'CAPE weighs as little as 19 pounds and features a versatile frame available in over 600 color combinations.





## ULTRA Lightweight Frame

Weighing only 19 lbs, the X'CAPE is one of the lightest pediatric folding wheelchairs available. This means energy that would have been spent on propulsion can be spent on playing, exploring and being a kid!

## Crossbrace with Patented XLOCK TECHNOLOGY

X'CAPE offers a level of performance and durability comparable to a rigid frame wheelchair. Its lockout feature and fewer moving parts reduce weight, frame play, and maintenance, ensuring more propulsion energy goes directly to the wheels. This allows your child to access all of their environments while experiencing less fatigue.

## TRUE Pediatric Design

The X'CAPE's setup – including seat widths, frame proportions, center-of-gravity settings, wheel sizes, and wheel angles (or camber) – can be fine-tuned to provide the best possible access to the wheels for each child. As they grow, 2" of width growth in the crossbrace and 3" of depth growth in the seat rail keep your child positioned properly for maximum independence.



ZIPPIE BELIEVES  
PROPER SEATING  
AND POSITIONING  
ARE KEY TO  
HEALTHY GROWTH  
AND DEVELOPMENT

## WC-19 Transit Approved

Approved for use as an occupant's seat in a motor vehicle when used with an approved Wheelchair Tiedown and Occupant Restraint System in accordance with SAE (J2249).



### One-step Push-to-Lock

One-step push-to-lock tab locks out the frame



### One-step Pull-to-Fold

One-step pull-to-fold strap collapses the frame effortlessly.



### Folding Push Handle


Tuck away behind the back upholstery or backrest while your child is self-propelling, or raise to your height so you can give an extra push.





### VersaRail™ Seat Rail

Pediatric wheelchairs with short seat rails often don't provide enough space for multiple components – custom seating, armrests, wheel locks, side guards, etc. – leaving you in the tough spot of choosing what is most important. VersaRail's unique design accepts both bolt-on mounts and mounting clamps, and allows for shared placement of these components, so your child has access to all the options they need—no compromises.

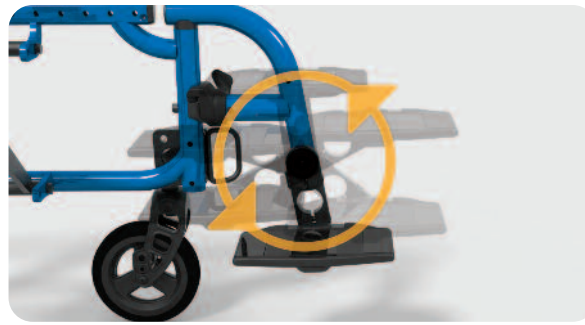

### Z-finity™ System

As your child grows and his/her needs change, the Z-finity System can be adjusted between .5" above the seat pan to 21" below the seat pan to provide proper support of the legs and feet at all times. Its rotating hardware allows infinite adjustment of the footrest's angle and placement to achieve optimal positioning of the feet, ankles, and knees.

EVERY DETAIL  
OF THE X'CAPE  
WAS DESIGNED FROM  
THE GROUND UP  
TO HELP CHILDREN  
TRANSITION TO  
THEIR NEXT LEVEL OF  
INDEPENDENCE



### JAY® Seating Options

VersaRail makes ordering simple by offering one seat rail system that is compatible with JAY Zip seating, JAY custom seating, and sling upholstery.



### Dynamic Back Option

Allows controlled movement of the backrest to help maintain proper user positioning and reduce wheelchair damage.



### Swing-In/Swing-Out & Fixed Hanger Options

Choose Swing-in/Swing-out hangers for aid in transfers and pulling up close to cabinets and counters. Choose Fixed hangers for the lightest weight configuration and an active look.



IMPORTANT NOTE: Print processes cannot accurately represent paint finishes. Actual colors may vary.



## Color Options

Choose from 24 standard frame colors and three special order frame colors. The crossbrace color is chosen separately at no extra charge to create over 600 custom color combinations!



## Spoke Guards



HCPCS Code: E1238

Transit Option: Yes\*

**Chair Weight:** 19 lbs  
8x8" fixed front frame, 20" lite spoke wheels. Excludes armrests, push handles & anti-tips. Actual weight may vary according to chair configuration & accessories.

**Weight Capacity**  
165 lbs (w/o seating)

**Overall Width**  
14.8" to 26.5"

**Seat Widths**  
8" to 16" (grows to 18")

**Seat Depths**  
8" to 18" (grows to 20")

**Tubing Diameter**  
Bottom Tube: 7/8"  
VersaRail: accommodates 1" clamps

**Back Angle Adjustment**  
Fore (-) 8° to  
Aft (+) 24°

**Front Seat-to-Floor Heights**  
13" to 19"

**Rear Seat-to-Floor Heights**  
11" to 18.5"

**Footrest/Frame Angle**  
80° (Zfinity™ may be adjusted to achieve a wide range of footrest angles)

**Back Types**  
Standard (adjustable) and Dynamic

**Back Heights**  
9" to 24"

**Back Post Options**  
Straight (no push handle),  
Integral Push Handle,  
Folding Stroller Handle

**Center of Gravity Adjust.**  
3" (-0.5" to 2.5")

**Vertical Axle Plate Adj.:** 4"  
**Rear Wheel Sizes**  
12", 16", 20", 22", 24"

**Armrest Options**  
T-Post w/Sideguard,  
Flip-up, Adjust. Locking  
Flip-up, Length Adjust.  
Locking Flip-up, Single  
Post Height Adjustable

**Caster Sizes**  
3", 4", 5", 6"

**True Fit or One-Time Free Growth Kit:** Yes

**Warranty:** Limited/Lifetime  
(Refer to Owner's Manual for details)

\*Conforms with ANSI/RESNA WC-19 standards. Refer to Quickie® & Zippie® WC-19 Transit Securement Supplement (MK-100072) for additional information. Supplement is available for download at [www.sunrisemedical.com](http://www.sunrisemedical.com).



©01.2015 Sunrise Medical (US) LLC  
MK-100117 Rev.B

800.333.4000 [www.SunriseMedical.com](http://www.SunriseMedical.com)

