



spex

Shaping Better Lives

Vigour
Cushion

Designed for the user who needs mild to moderate postural control.





30mm

Vigour Standard Contour



50mm

Vigour High Contour

Vigour Technology

Vigour cushion is designed to properly centralise and stabilise the pelvis, hips, and legs to maintain and improve comfort, function, and positioning in all other parts of the body.

Vigour reduces fatigue and pressure with the support and comfort that helps you go the distance. A combination of high resilient foam and gel-foam that is durable, flexible and comfortable.




The Vigour Cushion is an excellent choice for individuals who:

- Require pressure reduction and redistribution.
- Require pelvic stability.
- Need moderate control of posture and upper extremities.
- Need moderate control of abduction or adduction propensities.





Functional Benefits

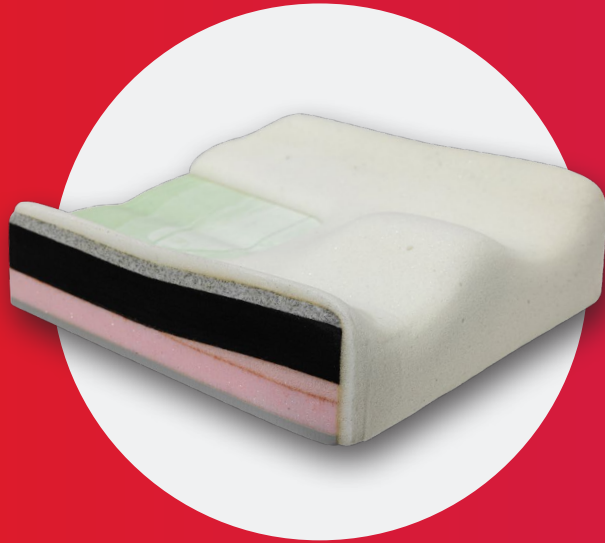

Gel-Foam Ischial Well

Evenly distributes pressure under the pelvis and in combination with the dual-layer cover, reduces shear forces in sensitive areas.



Built-in Thigh Ramp

The Vigour cushion features a strategically designed in-built thigh wedge at the front of the cushion, increasing user comfort and enhancing pelvic stability by optimising thigh loading.




Comes standard with dual-layer outer cover and incontinence cover

Dual-layer cover reduces shear forces, and breathable air-mesh on cushion sides counters moisture.

Accessories & Customisations

Colours Available:

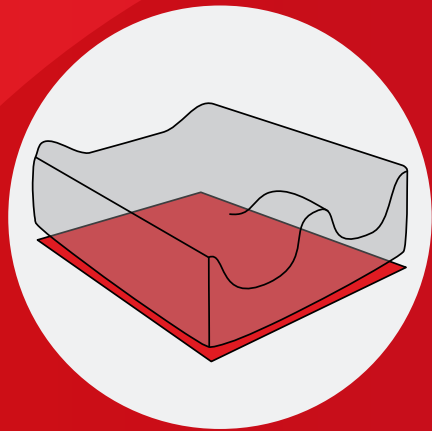


Cushion Rigidiser

Code: 1110-WWDD-200



**Spare Cushion
Outer Cover**

Code: 1209-2491-XXX

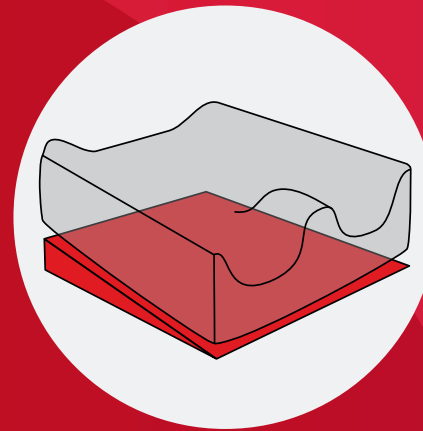


5° Wedge

Code: 1209-2724-300



**Phone
Accessory Bag**


Code: 1111-0480-000

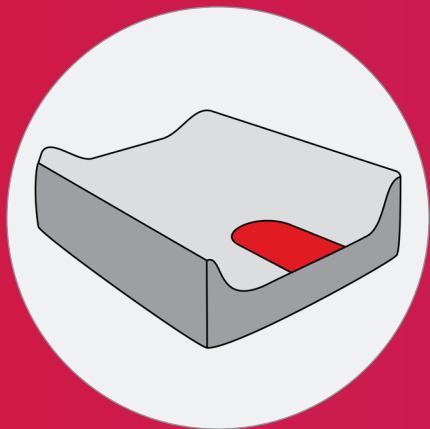

Remove Medial Thigh Contouring

Remove the Medial Thigh Support on the cushion so that a flat front surface is achieved without any contouring. This is ideal relief for excess tissue and a great option to facilitate easy transfers for users without natural adduction of the thighs.


Code: 1111-0335-000

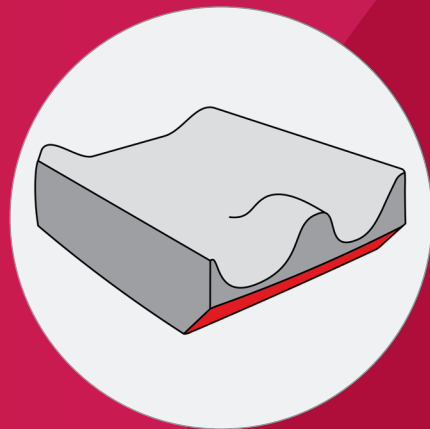

Front Surface Cut-back


Front surface cut-back, or a sloping front, can be achieved for users with excess tissue or to provide relief for users when transferring and/or foot propelling.

Code: 1111-0281-000

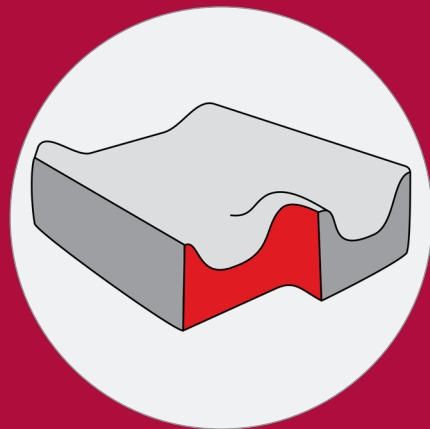

Leg Length Discrepancy

Code: 1111-0292-000

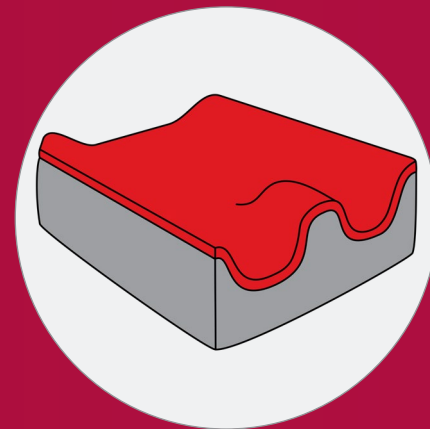

Additional Pressure Relieving Overlay

Code: 1111-0391-000

Supracor® Honeycomb Overlay


Code: 1111-0431-000

Refer to the Spex Build-Your-Own Cushion script form for a detailed list of customisations.

Specifications




CUSHION SIZE CHART

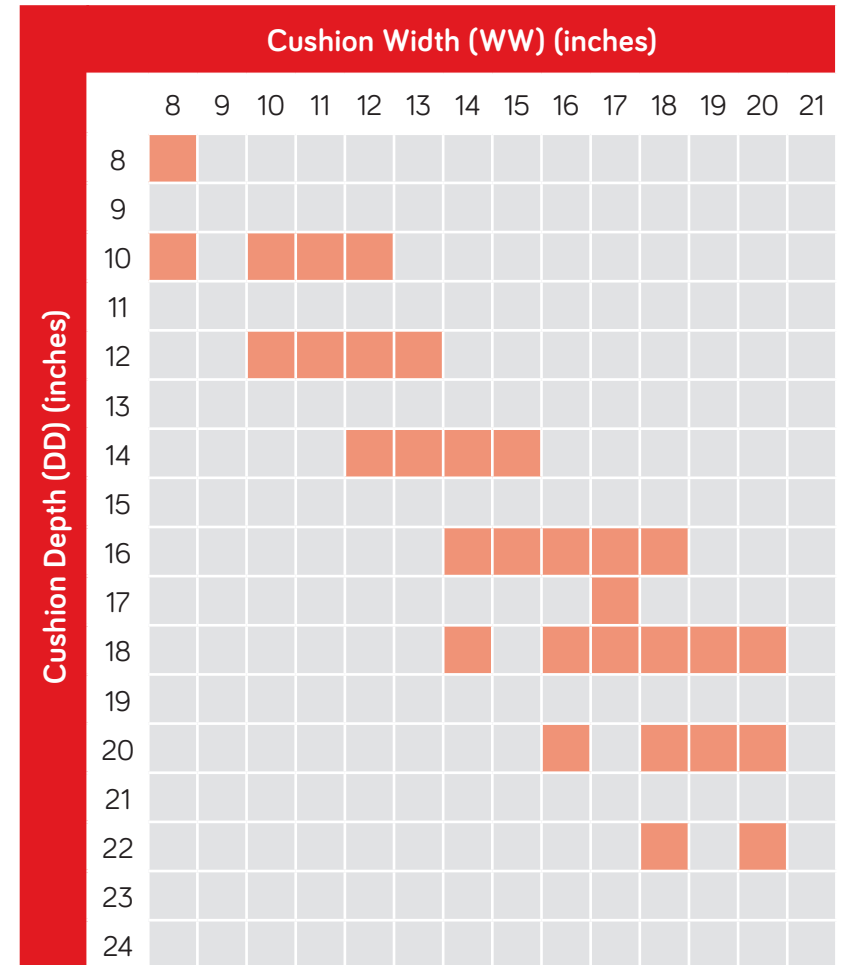

- Red squares are standard sizes.
- Non-standard or custom sizes are available on request.

CUSHION DETAILS

Width	Height	Capacity	Trough Height		Product Weight	
			Std	High	Std	High
8"-10"	2"	200lbs (90kg)	30mm	50mm	(16x16)	(16x16)
11"-15"	2.5"					
16"-17"	2.75"	350lbs (160kg)			1.4kg	1.5kg
18"-20"	3"					

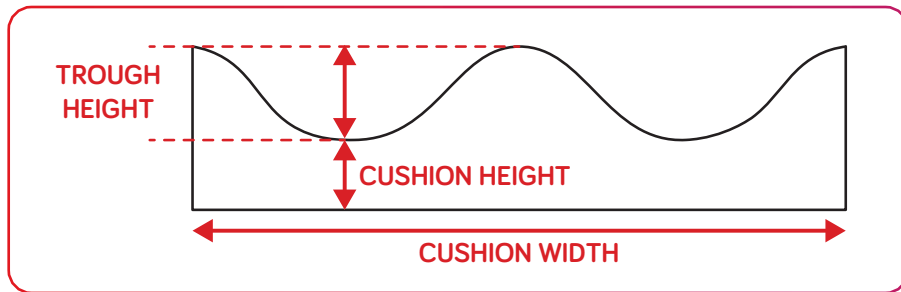

CUSHION ORDER CODES

Standard Contour	High Contour
1103-WWDD-OSC	1103-WWDD-OHC



PRODUCT PERFORMANCE

Pressure Care			3	
Positioning			3	
Lateral Stability			3	
Offloading			3	
Changing Needs	1			

- Note: This is a guide only and ratings are dependent on correct use of product.

spex

Manufactured by Spex Ltd
solutions@spexseating.com
spexseating.com

DISTRIBUTED BY:

Flyer: Spex Vigour Cushion_20240611

