

magic
Mobility™

Magic XT2

*Your Journey,
Your Way*



magic
XT2™

magic
XT2™



LIFE'S A *Journey*

The XT2 is designed to journey far and wide outdoors, to let you see and do more than ever before. Confident in the stability of a long wheelbase rear wheel drive with fully independent front and rear suspension. Powered by massive motors and long-range batteries. Manoeuvring with the ease and precision of servo-powered, joystick-controlled steering. Enabled by personalised seating options and equipped with features and accessories designed for your outdoor journey.

Feel free to journey further. It's *your* journey, so do it your way.

INTEGRATED LED LIGHTING

Don't get caught out in the dark! Light your way with confidence thanks to our integrated front and rear LED lighting system.



POWER THROUGH WITH EASE

1,700 Watts of power

Each of the XT2's rear wheels are driven by powerful 850 Watt gear in line motors, delivering plenty of power to zip around, up, or over and reach the places you want to go.

PATENTED SUPERIOR SUSPENSION

Independent front & rear suspension

Each wheel is mounted on its own pivoting arm that adjusts independently to follow the contours of your pathway. Keeping all four wheels firmly on the ground even while driving at higher speeds over uneven surfaces. Delivering greater stability, comfort and control for a feeling of effortless driving confidence.

ALL-DAY RANGE

115 Amp Hour batteries

High capacity batteries come standard on the XT2, providing increased running time on a single charge over our standard batteries. Experience the range you need to enjoy an extended outdoor journey knowing you'll safely make it home at the end of a long and limitless day.

*115Ah at C20 discharge rate.



PRECISE FORWARD & REVERSE CONTROL

Servo powered steering

The XT2 steers using a combination of differential wheel speeds and an industry-leading servo steering module. The chair's intuitive joystick steering delivers a feeling of reliability, responsiveness, and effortless control at any speed.

Because the world's not flat™

JOURNEY IN COMFORT & STYLE

Equip your chair for the journey ahead. Choose from a wide range of options to meet all your needs for comfort, support and personal style.

OPTIONS FOR YOUR EXTENDED JOURNEY

Equip your chair with all you need for a day-long journey! Choose accessories such as cup holders, storage boxes, luggage racks, accessory bags, phone mounts and chargers for convenience. And select the power seating options and features you need to be comfortable all day.



REHAB OR MPS SEATING

Choose from our rehab or MPS seating and backrests, and a wide range of other options including fold-forward or removable backrests for ease of transport, swing away joystick mounts and a wide range of headrests. All core seating options are standard across the Magic Mobility portfolio.



CHOOSE YOUR BASE COLOUR

Select from our three most popular colours. Our sleek, high-quality metallic finish adds wow factor, and gives depth and dimension to the XT range.



XT2 SPECIFICATIONS

| | | |
|---------------------------|--|--------------------------|
| USER WEIGHT LIMIT | | 136kg / 300lb |
| MAXIMUM SPEED | | 10km/h / 6.2mph |
| MAXIMUM SPEED (with lift) | | 5km/h / 3.1mph |
| LENGTH OF BASE | | 1105mm / 43.5" |
| WIDTH OF BASE | | 685mm / 27" |
| WHEELS | | Urban pneumatic 14" x 3" |
| SEAT TO FLOOR HEIGHT | | 480mm / 19" |
| GROUND CLEARANCE | | 125mm / 5" |
| TURNING RADIUS | | 1750mm / 69" |
| MOTORS | | 2 x 850W gear in line |
| CHAIR WEIGHT (minimum) | Base, batteries, rehab seating, centre post legrest, no power lift functions | 177kg / 390lb |
| (standard) | Base, batteries, rehab seating with power lift, tilt, recline and legrest | 212kg / 467lb |
| DRIVE ELECTRONICS | | R-net |
| BATTERY SIZE | | 115Ah (C20) |
| SEATING (width & depth) | | 305-610mm / 12-24" |

STANDARD POWER FUNCTIONS


Seat lift and tilt
300mm / 45°

95-180°
Legrest
elevation

Tilt 50°

95-160° Recline

Seat lift
12" / 300mm

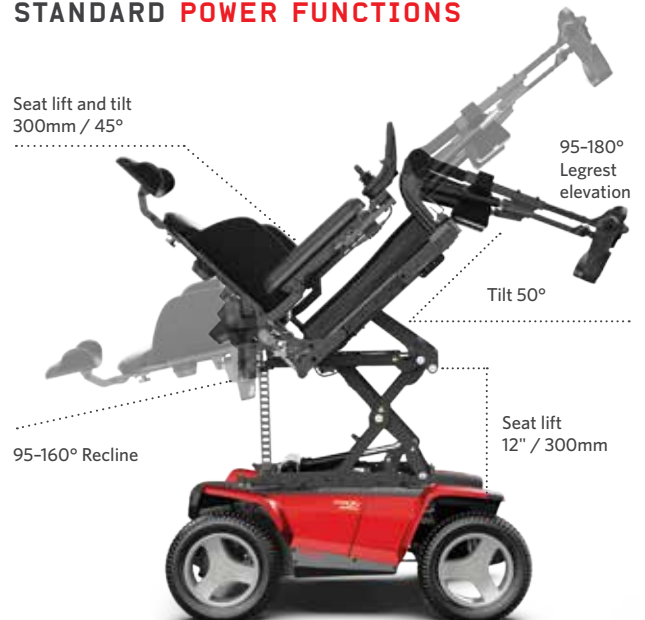

Because the world's not flat



Magic Mobility +61 3 8791 5600
sales@magicmobility.com.au
Explore with us magicmobility.com.au



Scan QR to learn more

Please remember that your wheelchair has electric motors and must NEVER be driven through water, into rivers, creeks and the sea. Check online for current version. The information and specifications contained in this brochure are correct and up-to-date at the time of printing and may be subject to change. ©Magic Mobility. Product of Australia.

